



RAVEN-SB SYSTEM SOFTWARE

RESIDENT ASSESSMENT VALIDATION AND ENTRY (RAVEN-SB) FOR SWING BED SYSTEM REFERENCE MANUAL

BETA VERSION 1.0

MAY 1, 2002

TABLE OF CONTENTS

OVERVIEW.....	1
<i>Introduction.....</i>	<i>1</i>
<i>Conventions.....</i>	<i>1</i>
QUICK START.....	3
<i>Installation and Start Up.....</i>	<i>3</i>
Hardware/Software	3
Printing RAVEN-SB Assessment Forms.....	3
<i>Install the RAVEN-SB System Application</i>	<i>4</i>
HELP.....	5
<i>Log In.....</i>	<i>6</i>
USER IDs and Passwords	6
Authorization Rights	7
SYSTEM SETUP.....	8
<i>Maintain Employee Database</i>	<i>8</i>
Add a New Employee Record	9
Edit an Employee Record.....	9
Save an Employee Record	9
Delete an Employee Record	9
Cancel an Employee Record.....	10
Authorization Rights	10
Allow Multiple Logins	10
Reset User Logon Flag.....	10
<i>Maintain Facility Database</i>	<i>11</i>
Facility Information.....	11
Medicare and State RUG Calculation Options	12
RUG Calculations	13
Agent Information	14
Facility/State Options	14
<i>Maintain Resident Database</i>	<i>14</i>
Add a New Resident Record	16
Edit a Resident Record.....	16
Save a Resident Record.....	17
Hide a Resident Record	17
View an Archived Resident Record.....	17
Delete a Resident Record.....	17
Cancel a Resident Record.....	17
<i>Enable Printing</i>	<i>17</i>
<i>Enable Backup</i>	<i>18</i>
<i>Windows Settings.....</i>	<i>18</i>
RAVEN-SB DATA ENTRY.....	19
<i>RAVEN-SB Management Screen</i>	<i>19</i>
Under the File Menu, the User Can:.....	20
Under the Options Menu, the User Can:	20
Under the Administration Menu, the User Can:	23

Under the Help Menu, the User Can:	29
Select/Add a Resident	29
Archive a Resident.....	30
Restore a Patient.....	30
Select/Add an Assessment.....	31
Special Features (RAVEN-SB Management Screen)	33
EXIT	33
RAVEN-SB Data Entry Window.....	34
Enter Data for a Variable	35
Special Features (RAVEN-SB Data Entry Screen)	37
Calculate/View RUGs (if Applicable)	39
SB_SUB_REQ.....	40
Assessment Validation after Exit.....	40
Assessment Status Window.....	41
RAVEN-SB DATA EXPORT.....	44
Exporting Assessments.....	44
Re-Exporting Assessments.....	45
RAVEN-SB DATA IMPORT.....	46
Selecting Import Options.....	46
Selecting to Show Warning Messages	47
Performing Error Checks on an Import File	47
Importing RAVEN-SB Data	47
RAVEN-SB HELP.....	49
General/About Help.....	49
Additional Help.....	49
REPORTS.....	50
Select a Report to View/Print	50
Assessment Data Entry Report.....	50
Report on an Assessment(s) Using the Resident Function.....	51
Report on an Assessment(s) Using the Status Function	52
Event Tracking Report	52
Report View Window.....	55
RAVEN-SB Assessment Forms.....	55
APPENDIX A.....	56
Repair/Compact Utility for Microsoft Access Database	56
Repair and Compact	56
GLOSSARY OF TERMS.....	58

OVERVIEW

INTRODUCTION

The Swing Bed Minimum Data Set (SB-MDS) was developed by the Centers for Medicare and Medicaid Services (CMS) to establish a standardized protocol for assessing the clinical characteristics of swing bed facility residents. This information will be used to better support the regulatory process and policy-making by providing survey agencies with a mechanism for objectively measuring and comparing facility performance and quality. In addition, it will provide researchers with a rich set of information to support the development of improved standards of care through the study of resident care outcomes.

CMS has sponsored the development of a data entry system known as the Resident Assessment Validation and Entry for Swing Bed System (RAVEN-SB). RAVEN-SB can be used by swing bed facilities for collecting SB-MDS data in a standardized format and electronically sending SB-MDS data records to the State System. Based on CMS's MedQuest tool, RAVEN-SB enables the user to electronically enter resident assessment data into SB-MDS standard record formats, calculate Resource Utilization Groups (RUGs), and export the resident assessment data records to the State System.

CONVENTIONS

Standard conventions have been adopted in the RAVEN-SB Reference Manual to facilitate locating and identifying desired information. These conventions and their definitions include the following:

Bolded Title Case Windows, screen titles, program functions, variable names, and dialog boxes (e.g., the **Resident Information** screen).

Title Case and Italicized Text box or list box descriptions (e.g., the *Select an Employee* list box).

Title Case Menu options (e.g., the General Help menu option).

<CAPITALS> Menus, buttons, icons, and function keys (e.g., the HELP menu, the <CANCEL> button, the <NOTES> icon, and the <T> key).

Bolded User entered data (e.g., Enter the password **SWB**).



Notes Symbol. This symbol indicates that the user should take notice of the information.



Question Mark Symbol. This symbol prompts a user to review the procedure about to be invoked. For example, if a user adds a new patient and then decides to cancel the addition, a message box with a question mark will appear asking the user to verify the request.



X Symbol. This symbol indicates that the requested procedure cannot be completed and displays a short message explaining why.



Exclamation Point Symbol. This warning symbol indicates that an action has been selected that may have undesirable results. The warning provides an explanation and asks for verification before continuing with the selected action.

QUICK START

INSTALLATION AND START UP

This section discusses the following topics:

- Hardware
- Software
- Printing
- Installation
- Help

HARDWARE/SOFTWARE

- IBM or IBM-compatible personal computer with Pentium 200 MHz or higher processor
- VGA or VGA-compatible display monitor with a minimum resolution of 800x600 pixels

NOTE:

RAVEN-SB was designed using Super VGA resolution (800x600). At a lower resolution (e.g., 640x480), all of the variables on a given screen may not be visible. Use the scroll bars to view the entire screen and perform data entry.

- Hard disk with a minimum of 30 megabytes of available space
- Sixteen megabytes of RAM (32 are recommended)
- Microsoft Windows 95, 98, NT, 2000, or XP
- Microsoft or Microsoft-compatible mouse or pointing device

PRINTING RAVEN-SB ASSESSMENT FORMS

Before printing the RAVEN-SB assessment forms, it will be necessary to activate the Adobe Acrobat Reader. If the software is not installed, the user can download the installation file from the Adobe web site (www.adobe.com/).

NOTE:

Instructions on how to activate the Adobe Acrobat Reader are provided in the System Setup section.

INSTALL THE RAVEN-SB SYSTEM APPLICATION

Detailed installation instructions are contained in README.TXT. Unzip the downloaded RAVEN-SB installation file in a temporary directory (e.g., C:\TEMP). This action will decompress all of the installation files as well as README.TXT. Then run the SETUP.EXE program to begin the installation and follow the steps to install the RAVEN-SB program.

Operating RAVEN-SB within Large Facilities

The intent of RAVEN-SB is to provide software to facilities to whom either: 1) could not afford to purchase software from a vendor; and/or 2) did not require a full-featured, client/server application that would be integrated with existing systems. As such, RAVEN-SB is designed as a stand-alone application and is not intended to be run in a multi-user environment. Consequently, RAVEN-SB utilizes Microsoft Access databases for storage. This places limits on the software in terms of performance and capacity.

However, RAVEN-SB will operate correctly in a network environment and allow entry and maintenance of a large number of assessments:

- The RAVEN-SB storage database can be located on a shared, network drive while the RAVEN-SB application files are located on the client PC. This capability allows multiple installations of the software to perform data entry against the same database. For more information, see the section below titled "Installing RAVEN-SB in a Network Environment".
- The Archive functionality is included in RAVEN-SB. This functionality alleviates the burden on the primary storage database and allows RAVEN-SB to operate in a high capacity environment. However, this functionality does not operate automatically and requires conscious effort on the part of the System Administrator. For more information, see the "Archive Setup" section.

RAVEN-SB should not be mistaken for a network application. It is still a stand-alone application, but contains features that make it possible to operate in a network environment under limited circumstances. Microsoft Access is not a full-featured RDBMS. Although a number of steps have been taken to address collisions in the software, collisions that will cause application errors can still occur. In addition, when the number of records in the storage database becomes very high, the performance and reliability of the system will degrade.

To ensure the best possible operation under a network environment, the following steps should be taken:

- Limit the number of collisions against the database
- Assign one user to perform all data entry for a particular resident
- Assign one user to perform all maintenance of the resident information including the adding and deleting of residents
- Assign one user to perform all maintenance of the employee information
- Assign one user to perform all maintenance of the facility information
- Assign one user to perform all imports and exports
- Assign one user to perform all archiving and restoring

- Do not perform maintenance of the resident information while other users are performing data entry
- Do not perform maintenance of employee information while other users are performing data entry
- Do not perform maintenance of facility information while other users are performing data entry
- Do not perform imports and exports while other users are performing data entry
- Do not perform archiving and restoring while other users are performing data entry
- Utilize the **Archive** functionality regularly to alleviate the burden on the storage database. In addition, be certain to compact the storage database (SWB.MDB) after several archives.

Installing RAVEN-SB in a Network Environment

There are several steps to follow when installing RAVEN-SB in a network environment. These steps include:

- For the first PC, install the RAVEN-SB software. Log in to the RAVEN-SB system and click on the ADMINISTRATION menu and select the Preferences option. Click on the Location tab and enter the network path where the SWBMDB and SWBTRACK.MDB databases will reside. Exit RAVEN-SB.
- Copy SWB.MDB and SWBTRACK.MDB to the network path that was specified in the previous step.
- Log in to RAVEN-SB and set up employee, facility, and resident information. If desired, specify an archive database. Also, go to the **RAVEN-SB Preferences** screen and select the desired options and paths. Exit RAVEN-SB and copy SWBDICT.MDB (in the installation directory) and SWB.INI (in the Windows SWBPRINT directory) to the network.
- Now install RAVEN-SB on all additional client PCs. Copy the SWBDICT.MDB from the network to the installation directory and copy the SWB.INI to the Windows SWBPRINT directory. This process will ensure that the settings established from the original installation are incorporated into each of the additional client installations.


HELP

The General Help provides RAVEN-SB system help regarding how to use and maintain the RAVEN-SB software. The About (RAVEN-SB) Help identifies software and dictionary version information. From the RAVEN-SB Management screen or any of the RAVEN-SB Data Entry screens, click on the HELP menu and select either the General Help or the About Help options.

NOTE:

Context sensitive help is available from any window by pressing the <F1> key. Individual variable help is also available by right mouse clicking anywhere within the variable's frame.

LOG IN

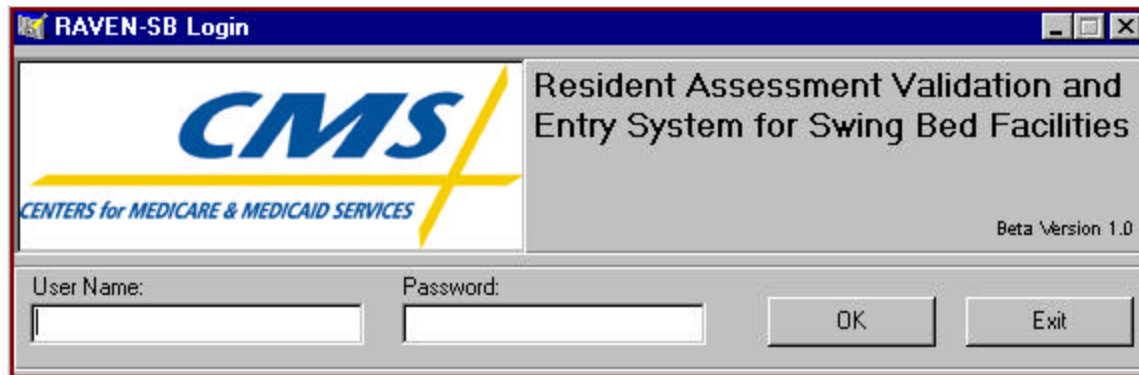

The user can log in to the RAVEN-SB module by entering the user identification and password in the *User Name* and *Password* text boxes, respectively.

USER IDs AND PASSWORDS

Use the following log in identifications and passwords for RAVEN-SB the first time:

- User ID **SWB** and password **SWB**

The user login identification and password can be added, modified, or deleted by using the **Maintain Employee Database** feature within the RAVEN-SB software. For information on Authorization Rights, see Authorization Rights under the System Setup section.

NOTE:

After installation, it is recommended that the System Administrator log in first, using **SWB** as both the User ID and Password. The first task that the System Administrator should complete is to establish a new User ID and Password that will be used as the System Administrator account. This task can be accomplished by using the **Maintain Employee Database** function described in Chapter 3, System Setup. The **SWB** account should then be deleted to ensure the security of the RAVEN-SB software.

AUTHORIZATION RIGHTS

The *Authorizations* radio buttons on the **Maintain Employee Database** screen -- Data Entry, Clinical Supervisor, and System Administrator -- determine the level of access for the currently selected employee. The following table displays the Authorization Rights for users by function.

<u>Item</u>	<u>Data Entry</u>	<u>Clinical Supervisor</u>	<u>System Administrator</u>	<u>View Only</u>
Maintain Employee Information Database	No	No	Add/Edit/Delete	No
Maintain Facility Information Database	Edit in Data Entry mode	Add/Edit	Add/Edit	No
Maintain Resident Information Database	Edit in Data Entry mode	Add/Edit	Add/Edit/Hide/Delete/Archive	No
Assessments	Add/Edit	Add/Edit	Add/Edit/Hide/Delete	View
Perform Data Entry for All RAVEN-SB Forms	Yes	Yes	Yes	No
Review Assessments Performed by Other Users	No	Yes	Yes	No
View Hidden Assessments	No	Yes	Yes	No
Correct Assessment	No	No	Yes	No
View User List	Yes	Yes	Yes	Yes
Calculate RUGs	Yes	Yes	Yes	No
Print Assessments	Yes	Yes	Yes	Yes
View/Print Reports	Assessment Data Entry Report	Assessment Data Entry Report	Assessment Data Entry Report Event Tracking Report	Assessment Data Entry Report
Set Tracking of Events	No	No	Yes	No
Set RAVEN-SB Preferences	No	No	Yes	No
Import Files	No	No	Yes	No
Export Files	No	No	Yes	No

SYSTEM SETUP

MAINTAIN EMPLOYEE DATABASE

The Employee Database provides centralized management of employee records that will be used continuously in the entry of assessment data and to customize the operation of the RAVEN-SB software. It is accessed from the main menu by selecting the ADMINISTRATION menu item and then choosing the Maintain Employee Database option.

Note that this screen is only accessible to users with System Administrator rights. After installation, the user (authorized as System Administrator) will want to define the expected users of the RAVEN-SB software. The **Employee Information** screen is used to collect employee identification and authorization information. This screen consists of two sections: Select an Employee and Enter Employee Information.

Employee Information

Select an employee:

	User ID	On	Login Time	Location	Phone	Resident	Assessment
1	SWB	*	16:11:17				

Enter Employee Information

User ID: Password: Title:

First Name: Middle Initial: Last Name:

Authorizations: ☐ Data Entry ☐ Clinical Supervisor ☐ Systems Administrator ☐ View Only

☐ Allow multiple logins

The *Select an Employee* section consists of: User ID; On * (indicates the employee is currently logged in); Login Time; Location; Telephone; current Resident and Assessment that the employee is working on (if the user is logged in), and the Status of that assessment for all employees currently defined in the database.

When the user clicks on a specific User ID, the D is highlighted and the corresponding employee data is displayed in the *Enter Employee Information* section. The *Enter Employee Information* section consists of:

User ID; Password; Title; First and Last Name and Middle Initial; Authorizations; Allow Multiple Logins, and Reset User Logon Flag.

The user may:

- Add a new employee record
- Edit an existing employee record
- Save all changes made to an employee record
- Delete an employee record
- Cancel all edits to an employee record
- Assign Authorization rights
- Allow user to login multiple times
- Reset user logon flag
- Reset user logon status

ADD A NEW EMPLOYEE RECORD

Click on the <NEW> button to add a new employee record. The User ID and Password fields must be completed. Upon completion of the Employee Information section, click on the <SAVE> button to update the employee database. Saving will also refresh the *Employee Selection* list box to display the new user.

EDIT AN EMPLOYEE RECORD

Click on the field or checkbox of interest to edit the data. Clicking the <SAVE> button will save all changes made to the employee record to the database, while clicking the <CANCEL> button will abandon all changes made to the employee record since the last save.

SAVE AN EMPLOYEE RECORD

Click on the <SAVE> button to save all information entered for a new or existing employee record.

DELETE AN EMPLOYEE RECORD

Click on the <DELETE> button to remove the currently selected employee record from the database.

NOTE:

Since the System Administrator has access to all functions of the RAVEN-SB software, the user is prevented from deleting the last System Administrator record from the database or changing the rights (from System Administrator to Data Entry) for the last System Administrator record.

CANCEL AN EMPLOYEE RECORD

Click on the <CANCEL> button to abandon all changes made to the employee record since the last save.

AUTHORIZATION RIGHTS

The *Authorizations* radio buttons on the **Employee Information** screen -- Data Entry, System Administrator, and Clinical Supervisor -- determine the level of access for the currently selected employee.

- **Data Entry Rights.** Users are able to edit Resident/Facility information (in **Data Entry** mode) add/edit assessments, perform data entry for all RAVEN-SB forms, calculate RUGs, print assessments, and view/print reports.
- **Clinical Supervisor Rights.** Users are able to add/edit Facility and Resident database information, add/edit assessments, review assessments performed by other users, view users, view hidden assessments, perform data entry for all RAVEN-SB forms, calculate RUGs, print assessments, view/print reports, and track resident/facility/assessment/employee information.
- **System Administrator Rights.** Users are able to add/edit/delete Employee database information, add/edit Facility database information, add/edit/hide/delete/archive Resident database information, add/edit/hide/delete assessments, create corrections, review assessments performed by other users, view users, view hidden assessments, perform data entry for all RAVEN-SB forms, calculate RUGs, print assessments, view/print reports, track resident/facility/assessment/employee information, set RAVEN-SB preferences, and import and export files.
- **View Only Rights.** Users are able to view what has been entered in data entry. Users are able to view/print reports, print assessment(s) from file, and view assessment data.

ALLOW MULTIPLE LOGINS

The *Allow Multiple Logins* check box allows the user to login using the same User ID from multiple installations in a network environment. For stand-alone installations, it is recommended that this option be checked.

RESET USER LOGON FLAG

This feature allows a System Administrator to reset the selected Employee's login status. After selecting an Employee and clicking on the <RESET USER LOGON FLAG> button, the Employee's status will be reset from "On" (i.e., an asterisk indicates that the user is logged in) to "Off", (i.e., the asterisk is removed). This function is useful when the RAVEN-SB software has not been exited properly such as after a power outage.

MAINTAIN FACILITY DATABASE

The Facility Database collects facility identification and system preferences and options. It is accessed from the main menu by selecting the ADMINISTRATION menu item and then choosing the Maintain Facility Database option.

The **Facility Information** screen collects facility, agent, and system information used to customize the operation of the RAVEN-SB software. It is recommended, after installation, a user with System Administrator rights complete this information. This screen is divided into five sections: 1) Facility Information; 2) Medicare and State RUG Options; 3) State RUG Options; 4) Agent Information; and 5) Facility/State Options.

The screenshot shows the 'Facility Information' window with a blue title bar. It contains five main sections:

- Enter Facility Information:** Includes fields for 9a. Facility Medicaid Provider Number, 9b. Facility Medicare Provider Number, Facility ID, Facility Name, Facility Address, Zip Code, City, State, Facility Contact Person Name, and Facility Contact Person Phone Number and Ext.
- Enter Agent Information:** Includes fields for Agent Tax ID, Agent Name, Agent Address, Zip Code, City, State, Agent Contact Person Name, and Agent Contact Person Phone Number and Ext.
- Medicare RUG Options:** Includes a Calculation Method section with radio buttons for Medicare (selected) and None, and a Case Mix Index dropdown menu set to 'A01 - Medicare PPS Rural'.
- State RUG Options:** Includes a Calculation Method section with radio buttons for Medicare (selected), Index, Hierarchical, and None, and a Case Mix Index dropdown menu set to 'A02 - Medicare PPS Urban'.
- Enter Facility/State Options:** Includes a checkbox for '100% Medicare/Medicaid certified facility'.

At the bottom, there is an 'Edit CMI Sets' button and 'Save' and 'Cancel' buttons on the right.

FACILITY INFORMATION

The Facility Information section enables the user to specify information that can uniquely identify the facility such as facility name, address, contact person, and Medicare, Medicaid, and State-specified provider numbers. Fields shown in **bold** are required to complete and export RAVEN-SB assessments.

NOTE:

When the user enters a valid zip code for the Facility Address, the correct city and state will be added. The user can still change the city and state if needed.

MEDICARE AND STATE RUG CALCULATION OPTIONS

RUG calculations may be performed for both Medicare and State purposes. Note that the user must specify at least one type of RUG calculation other than None (Medicare, Hierarchical, or Index) for either Medicare or State purposes when completing the **Facility Information** screen.

The user has the option to specify the RUG calculation method (Medicare, Hierarchical, Index, or None, if RUGs will not be calculated) for both the Medicare (43MDCR) and State (43STATE) variables in SB-MDS assessments. These methods are described further in the RUG calculation pseudo-code documentation available on the CMS SB-MDS web site.

For each RUG calculation to be performed, the user must also select a CMI Index set. Select a CMI Index set from the Case Mix Index drop down menu list on the **Facility Information** screen. Three standard sets are available for selection in this release:

- Medicare PPS Rural (A01)
- Medicare PPS Urban (A02)
- Research 44 Group (B02)

Note that these index sets are read-only; they cannot be deleted. In addition, the individual index values for the standard sets cannot be modified. Click on the <VIEW CMI VALUES> button to view the index values for the highlighted CMI set. After viewing the group index values, click on the <CLOSE> button.

Edit CMI Set for RUG Calculations

Medicare PPS Rural--07/01/98
Medicare PPS Urban--07/01/98
Research 44 group--Nurs. only
Test--

Add Delete View CMI Values

Description Edit Save

Code Date

Exit

In RAVEN-SB, new CMI sets are created from one of the three standard sets; this feature is useful if the new set only requires changing one or two of the index values from the standard set. To create a new CMI Index set, click on the <EDIT CMI SETS> button from the **Facility Information** screen. Highlight the appropriate set and click on the <ADD> button. This action will create a working copy of the standard index set that the user can name and modify as needed. Use the <EDIT> button to modify the description for the new index set and the <SAVE> button to store it in the database.

NOTE:

After a new index set has been added, the context sensitive <VIEW CMI VALUES> button will change to <EDIT CMI VALUES>. If a read-only set of indices is highlighted, the button will read <VIEW CMI VALUES>. If the user has created the highlighted index, the button will read <EDIT CMI VALUES>.

A user-created index set can be deleted by highlighting the index set and clicking the <DELETE> button. As noted above, the standard CMI Index sets are read-only and cannot be deleted.

	GROUP	CMIINDEX
1	RUC	44.000
2	RUB	43.000
3	RUA	42.000
4	RVC	41.000
5	RVB	40.000
6	RVA	38.000
7	RHC	39.000
8	RHB	36.000
9	RHA	33.000
10	RMC	37.000
11	RMB	34.000
12	RMA	32.000
13	RLB	31.000
14	RLA	26.000
15	SE3	35.000
16	SE2	30.000
17	SE1	29.000

Click on the <EDIT CMI VALUES> button to change the CMI values as needed. The user can save the changes by clicking on the <SAVE> button or exit the **Edit CMI Set** screen by clicking on the <CANCEL> button.

RUG CALCULATIONS

RUG calculations are automatically performed for an assessment when the data entry for the assessment has been completed and its status is made "Export Ready". RUGs can be calculated for all applicable assessments. The RUGs can be recalculated at any time during data entry. Click on the RUGs menu item in data entry to recalculate the RUGs. The **Calculate RUGs** screen will be displayed. Note, however, that the RUG values are not saved until the assessment is locked.

RUGs can also be calculated for an applicable assessment from the **RAVEN-SB Management** screen using the <RUGS> button. If the assessment is "Export Ready" or "Exported", the RUGs have already been calculated and will merely be displayed in the **View RUGs** screen. RUGs will be recalculated and displayed for applicable assessments with any other status. If RUGs are not applicable for an assessment (e.g., a

discharge or reentry assessment), then the user is notified that RUGs cannot be calculated for that assessment type.

AGENT INFORMATION

The Agent Information section collects unique agent data such as name, address, tax identification number (TIN), and telephone number.

NOTE:

The Facility Information and Agent Information data are included as header information for each export from the RAVEN-SB software. Information about RAVEN-SB export capabilities is provided in the Export section.

When the user enters a valid zip code for the Agent Address, the correct city and state will be added. Users can change the city and state.

FACILITY/STATE OPTIONS

If the facility contains only Medicare—and/or Medicaid—certified beds, then check this option. This action sets the SB_SUB_REQ to 1 by default for newly created assessments (not existing assessments). Note that the user still has the option to change SB_SUB_REQ within data entry, if desired.

MAINTAIN RESIDENT DATABASE

The Resident Database provides centralized management of resident information records that will be used continuously in the entry of assessment data and also to customize the operation of the RAVEN-SB software. It is accessed from the main menu by selecting the ADMINISTRATION menu item and then choosing the Maintain Resident Database option.

The **Resident Information** screen collects resident information that is reused in RAVEN-SB assessments. It is recommended that a user with System Administrator rights complete this screen prior to creating assessments for the residents. This screen consists of two sections: the Resident Selection section and the Resident Information section.

The Resident Selection section on this screen is similar to the Resident Selection section on the **RAVEN-SB Management** screen. The only option available in the Resident Selection section on the **RAVEN-SB Management** screen is the capability to add a resident. Any maintenance (e.g., edit/hide/delete) for a resident record must be completed through the Maintain Resident Database menu option.

Resident Information

Select a Resident

Search for... Last Name = [] Search Filter for... Unit = [] Filter

Last Name	First Name	MI	SSN	Unit	Medicare No.
[Empty Table]					

Browse Next 1000

☐ View Archived Residents Also

☐ View Hidden Residents Also

1. Resident Name

a. (First) [] b. (M.I.) [] c. (Last) [] d. (Jr./Sr.) [] ☐ Hidden

2. Gender

☐ 1. Male ☐ 2. Female

3. Birthdate []

4. Marital Status

☐ 1. Never Married ☐ 2. Married ☐ 3. Widowed ☐ 4. Separated ☐ 5. Divorced

5. Race/Ethnicity: (Check all that apply)

☐ a. American Indian/Alaskan Native ☐ b. Asian ☐ c. Black or African American ☐ d. Hispanic or Latino ☐ e. Native Hawaiian or other Pacific Islander ☐ f. White

6. Zip Code []

7. Resident SSN and Medicare Numbers

a. Social Security Number [] b. Medicare or Railroad Insurance Number []

8. Resident Medicaid Number [] [Enter '+' if pending or 'N' if not a Medicaid recipient in first digit followed by blanks]

State Resident ID Code []

Facility Specific Resident ID []

Unit []

Delete New Save Cancel

This section allows the user to select, search, or filter for a resident in the same manner as the **RAVEN-SB Management** screen. Residents are displayed in alphabetical order. Use the scroll bars on the *Resident* list box to view the available residents and click on the preferred row to select a resident. If there are more than 1,000 residents from which to choose, click on the <BROWSE NEXT 1000> button to display the next 1,000 residents.

The user may search the Resident database by selecting the Last Name, Social Security Number (SSN), Medicare Number, Medicaid Number, State Resident ID, or Facility Resident ID of a resident in the *Search for...* drop-down list box.

Select a Resident

Search for...

Last Name

Last Name

SSN

Medicare Number

Medicaid Number

State Resident ID

Facility Resident ID

After selecting an item from the *Search for...* drop-down list box, click on the <SEARCH> button to begin the search. The closest match will be the first item displayed in the *Resident* list box. To remove the search, clear the = text box and click on the <SEARCH> button. The residents will then be displayed alphabetically for selection within the *Resident* list box.

The user may want to reduce the number of residents available from which to choose by filtering on Unit. This filtering is accomplished by entering the desired Unit in the *Filter for Unit...* text box. Click on the <FILTER> button to begin the filter. An alphabetical list of residents by unit will be displayed in the *Resident* list box. If the Unit list is extensive, the user can search this list by Last Name, Social Security Number (SSN), Medicare Number, Medicaid Number, State Resident ID, or Facility Resident ID. To remove the filter, clear the Unit field and click on the <FILTER> button. An alphabetical list of residents will then be displayed for selection within the *Resident* list box.

The user can also sort by a field, (e.g., Last Name, First Name, MI, SSN, Unit, Medicare Number, Medicaid Number, State Resident ID, or Facility Resident ID) by clicking on the column header. Clicking once on the column header will sort the list in ascending order. Clicking on the column header again will sort the list in descending order.

Once a resident is selected, the corresponding data appears in the *Resident Information* text box. This section consists of data entry fields for demographic and identification information such as gender and race/ethnicity. The information collected is utilized to prefill the corresponding fields on the RAVEN-SB assessment forms for the selected resident.

The user may:

- Add a new resident record
- Edit an existing resident record
- Save all changes made to a resident record
- Hide a resident record
- View an archived resident record
- Delete a resident record
- Cancel all edits to a resident record

ADD A NEW RESIDENT RECORD

Click on the <NEW> button and begin entering data in the Resident Information section. The user must complete Last Name, Gender, Birthdate, and Race/Ethnicity before the record can be saved. The remaining fields are optional. Click on the <SAVE> button or the <CANCEL> button when data entry is complete. The user will be prompted with a message to "Save new resident currently displayed"? Click on the <YES> button to save the changes or the <NO> button to cancel the changes.

EDIT A RESIDENT RECORD

Click on the field or radio button of interest to edit the data. Clicking on the <SAVE> button will save all changes made to a resident record to the database, while clicking on the <CANCEL> button will abandon all changes made to the resident record since the last save.

SAVE A RESIDENT RECORD

Click on the <SAVE> button to save all information entered for a new or existing resident record.

HIDE A RESIDENT RECORD

Check the *View Hidden Residents Also* text box, select the resident to hide, and then check the *Hidden* text box. The resident is no longer visible in the *Resident* list box. Note that if a resident is marked as hidden, the user will not have access to the assessments for that resident. To view the hidden residents, check the *View Hidden Residents Also* check box. The user has the capability to redisplay a resident by unchecking the *Hidden* check box.

VIEW AN ARCHIVED RESIDENT RECORD

Check the *View Archived Residents Also* text box to view the archived residents. Until a resident is restored, the only activity that can be performed on an archived resident record is to search the list and view that resident record.

DELETE A RESIDENT RECORD

Click on the <DELETE> button to delete the currently selected resident. This action will remove the data record for the resident from the database. Note that a resident record can only be deleted if no assessments are associated with the resident.

CANCEL A RESIDENT RECORD

Click on the <CANCEL> button to abandon all changes made to the resident record since the last save.

ENABLE PRINTING

Before accessing the **Print** function for the RAVEN-SB assessment forms, it will be necessary to activate the Adobe Acrobat Reader. If this software is installed, then a System Administrator can click on the ADMINISTRATION menu item, select the Preferences option, and click on the Print tab. The System Administrator can then enter the file path/name into the *Adobe Acrobat Reader File Path/Name* text box (i.e., C:\PROGRAMS FILES\ADOBE\ACROBAT5\READER\ACRORD32.EXE). Completing the path for the Adobe Acrobat Reader from the PREFERENCES menu allows users to access the **Print** function for the RAVEN-SB assessment forms. Leaving the path blank prevents users from having access to the **Print** function for the RAVEN-SB assessment forms.

NOTE:

If the Adobe Acrobat Reader is not installed, then the user can download the installation file from the Adobe web site (www.adobe.com/).

ENABLE BACKUP

RAVEN-SB provides the option to make a backup copy of the data collection databases (e.g., SWB.MDB, SWBTRACK.MDB, EXPORT.MDB, RUGS.MDB, and SWBARCH.MDB each time the user exits RAVEN-SB normally). If these databases are corrupted due to abnormal termination of RAVEN-SB (e.g., power failure, system lockup), the System Administrator may copy the backup databases into the installation directory and only the last RAVEN-SB session will be lost.

To enable the backup, the System Administrator may click on the ADMINISTRATION menu item, select the Preferences option, and click on the Backup tab. Then the user may specify a directory location for the backup files. Note that this directory cannot be the directory where the application and databases currently reside. Also, if a path is not provided, the backup option is disabled.

NOTE:

The RAVEN-SB **Backup** capability is no substitute for a comprehensive backup strategy. It is **HIGHLY RECOMMENDED** that RAVEN-SB users institute a regular procedure of backing up the databases to a reliable media (e.g., floppies, CD-ROM, tape). For additional suggestions regarding backup, read the RAVEN-SB Help Desk Note on Backup (BACKUP.PDF).

Three compression options are available for the backup file (S). If "No Compression" is selected then the files are simply copied to the backup location. If "Enable Compression of Databases to ZIP Format" is selected then the backup files are compressed to SWBMDB.ZIP in the backup location. If "Enable Compression of Databases to Self-extracting EXE" is selected then the backup files are compressed to SWBMDB.EXE in the backup location. If an existing compressed backup file is present in the backup location then it is copied to SWBBACK.ZIP (or SWBBACK.EXE if the EXE option is selected) in the application directory.

WINDOWS SETTINGS

In order to ensure the reliable operation of the RAVEN-SB software, the following two Windows properties **must** be set accordingly:

- The Windows short date style must be set to display a four-digit year (e.g., 1999) instead of a two-digit year (e.g., 99).

For Windows 95, 98, NT, 2000, and XP click on the <START> button, select the Settings option and then select the Control Panel option to display the **Control Panel** screen. Next, click on the <REGIONAL SETTINGS> icon to display the **Regional Settings Properties** screen. Then, click on the Date tab to review the short date style properties. Verify that the short date style includes a four-digit year. A style that includes "YYYY" will satisfy this condition. The recommended setting is "MM/DD/YYYY".

- For Windows 95, 98, NT, 2000, and XP the Windows Font Size property must be set to "Small Fonts" and not "Large Fonts". This property can be found under the Settings tab of desktop properties.

RAVEN-SB DATA ENTRY

RAVEN-SB MANAGEMENT SCREEN

Once RAVEN-SB has been loaded, a number of functions may be performed. These functions are all activated from icons or drop-down menu items available on the screen.

RAVEN-SB Management

File Options Administration Help

Search for... Last Name = [] Search Filter for... Unit = [] Filter

Last Name	First Name	MI	SSN	Unit	Medicare Number
-----------	------------	----	-----	------	-----------------

Browse Next 1000 Archive Resident Add Resident

None Correct Assessment Add Assessment

Date	Type	#	Status	User	Open
------	------	---	--------	------	------

☐ View Only (Data Entry Not Allowed) ☐ View Hidden Assessments Also

Print RUGs Enter Data Exit

The **RAVEN-SB Management** screen consists of the following components:

- **Title Bar.** Displays the title **RAVEN-SB Management**.
- **Menu Bar.** Displays the list of functions that can be performed on the **RAVEN-SB Management** window. Each of these functions contains additional options. These options are displayed by using the mouse.
- **Button Bar.** Provides quick access to commonly used menu items. The functions available to all users include: **Log Off/On**; **User List**; **View/Print Reports**; **Print Assessment(s)**; and **Help**. An additional function, **Maintain Facility** and **Resident Database** information, is available to Clinical Supervisors and System Administrators. **Import**, **Export**, **Hide/Display Assessment**, **Hide All Assessments (Current**

Resident), **Delete Assessment**, and **Maintain Employee Database** information are available only to a System Administrator.

NOTE:

Additional features available to the System Administrator, but not displayed on the Button Bar, include: **Change SB_SUB_REQ**; **Archive Setup**; **Tracking Setup**; and **Preferences**. All of these features can be accessed from the ADMINISTRATION menu item.

- **Select a Resident.** Enables the user to select a resident to access his/her assessment history.
- **Add a Resident.** Allows the user (Clinical Supervisor and System Administrator) to add a resident.
- **Select Assessment for Resident.** Allows the user to view the history of previous assessments and their status for the currently selected resident. For additional information on the status of assessments, see the Assessment Status section below.
- **Add an Assessment.** Allows the user to add an assessment to a selected resident.
- **Correct an Assessment.** Allows the user to correct an assessment that has already been exported and accepted by the State System.

UNDER THE FILE MENU, THE USER CAN:

Open Archive

This option allows the System Administrator to open the archive directory and view the archived residents. To view an archive database, select the Open Archive option from the FILE menu on the **RAVEN-SB Management** screen. By opening the archive database, the user also has the option to restore a resident. See the Restore a Resident section for a complete explanation of the **Restore** function.

Exit

This option allows the user to exit the RAVEN-SB software. To exit RAVEN-SB, select the Exit option from the FILE menu on the **RAVEN-SB Management** screen or click on the <EXIT> button.

UNDER THE OPTIONS MENU, THE USER CAN:

Log Off/On

This option allows the user to return to the **RAVEN-SB Login** screen. From here, the user can log in as a different user or click on the <CANCEL> button and return to the **RAVEN-SB Management** screen.

User List

This option allows a user to view others who are logged on and the resident and assessment on which they are working. Select the User List option from the OPTIONS menu to display the **User List** screen. The **User List** screen contains the same information as the *Select an Employee* section of the **Employee Information** screen. This consists of: User ID; On (* indicates the employee is currently logged in); Login Time; Location; Telephone; name of current Resident and Assessment the employee is working on, and the Status of the assessment number.

User List							
	User ID	On	Login Time	Location	Phone	Resident	Assessment
1	SWB	*	12:44:54				

When using RAVEN-SB across a network environment, only one user can access a particular assessment at a time. The User List option displays any users who are currently logged on to the system and the assessment that user has open for data entry. Users can also see whether an assessment is in use by an asterisk (*) in the *Open* column for each assessment on the **RAVEN-SB Management** screen.

GERSHWIN						
			Correct Assessment		Add Assessment	
	Date	Type	#	Status	User	Open
1	04/10/2002	Non-PPS Assessment ()	00	New	SWB	*
2	04/10/2002	Discharge (Rtn not anticip..)	00	New	SWB	*
3	04/09/2002	PPS Assessment (5 day)	00	In Use	SWB	

NOTE:

An assessment is considered open while data entry is being performed. This situation only occurs when an assessment is accessed in a network installation of RAVEN-SB. The *Open* box will be empty when a user is working in a stand-alone configuration.

View/Print Reports

This feature is discussed in the Reports section.

Print Assessment(s) from File

This feature is discussed in the Reports section.

Import

This feature is discussed in the Import section.

Export

This feature is discussed in the Export section.

Hide/Display Assessment

This option gives the System Administrator the capability to specify whether an exported, modified, or inactivated assessment for the current resident can be viewed in the *Resident Assessment* list box. To hide an exported, modified, or inactivated assessment, highlight the desired assessment and select the Hide/Display Assessment option from the OPTIONS menu or click on the <HIDE/DISPLAY ASSESSMENT> button. The selected assessment will be removed from the list and given the status of "Hidden". Hidden assessments are only visible when "View Hidden Assessments Also" is checked. This feature is automatically enabled when an assessment is hidden.

To display a hidden assessment, click on the *View Hidden Assessments Also* check box to display the hidden assessment(s). Highlight an assessment and select the Hide/Display Assessment option from the OPTIONS menu or click on the <HIDE/DISPLAY ASSESSMENT> button. When the *View Hidden Assessments Also* check box is unchecked, the assessment will remain in the list. The assessment's status will be reset to the status it had prior to being hidden.

NOTE:

Although a Clinical Supervisor cannot hide/unhide an assessment, he/she can view hidden assessments by checking the *View Hidden Assessments Also* check box.

Hide All Assessments (Current Resident)

This option gives the System Administrator the capability to remove all exported assessments for the current resident from view in the *Resident Assessment* list box. Highlight a resident and select the Hide All Assessments (Current Resident) option from the OPTIONS menu or click on the <HIDE ALL ASSESSMENTS (CURRENT RESIDENT)> button. All assessments marked exported will be hidden from view and given the status of "Hidden".

To display a hidden assessment, click on the *View Hidden Assessments Also* check box to display the hidden assessment(s). Highlight an assessment and select the Hide/Display Assessment option from the OPTIONS menu or click on the <HIDE/DISPLAY ASSESSMENT> button. When the *View Hidden Assessments Also* check box is unchecked, the assessment will remain in the list. The assessment's status will reset to the status it had prior to being hidden.

NOTE:

Although a Clinical Supervisor cannot hide/unhide an assessment, he/she can view hidden assessments by checking the *View Hidden Assessments Also* check box.

Delete Assessment

This option gives the System Administrator the capability to delete a specific assessment for the current resident. Note that the assessment is deleted not only from the *Resident Assessment* list box but also from the database. Once an assessment has been deleted, it cannot be retrieved.

To delete an assessment, highlight an assessment and select the Delete Assessment option from the OPTIONS menu or click on the <DELETE ASSESSMENT> button to delete the assessment.

NOTE:

Assessments with a status of "Exported" or "Hidden" cannot be deleted. Once an assessment has been deleted, it cannot be retrieved.

UNDER THE ADMINISTRATION MENU, THE USER CAN:

The maintenance functions for the Employee, Facility, and Resident databases are discussed in the System Setup section. Users may also set up Archive and Tracking functions as well as specify RAVENSB preferences. Furthermore, an exported assessment requiring a manual change to SB_SUB_REQ can be modified using the Change SB_SUB_REQ feature.

Change SB_SUB_REQ

The Change SB_SUB_REQ function is used to modify the SB_SUB_REQ for an exported assessment that has been submitted to and accepted by the State System. Since the SB_SUB_REQ for this assessment cannot be changed on the State System by a modification request, it must be performed manually by State System personnel. Therefore, in order for the RAVEN-SB database to match the State System for this assessment, highlight the assessment and the Change SB_SUB_REQ function will be enabled since the assessment has Exported status. Select this option and the SB_SUB_REQ screen will appear. The user may then make the appropriate change. After clicking the <SAVE> button the new SB_SUB_REQ value will be stored and the assessment will maintain its Exported status.

Archive Setup

The **Archive Setup** function allows the user to define an archive database for a resident's assessments. By default, when the RAVEN-SB software is installed, the SWBARCH.MDB archive database is created.

The System Administrator has the capability to add an archive database. After selecting the ADMINISTRATION menu item, select the Archive Setup option to display the **Archive Information** screen.

Archive Information

Select an archive database:

	Name	Active	Location	Description
1	SWBARCH	*	C:\SWB\	The RAVEN-SB default

Enter Archive Information:

Database Name: .MDB

Location (Path):

Description:

☒ Make this database the active archive database.

The **Archive Information** screen displays two separate sections. The *Select an Archive Database* section displays the name of the archive database, a flag indicating whether that archive database is the active archive database, the location of the archive database, and a description. To add a new archive database, click on the <NEW> button and enter the following information in the *Enter Archive Information* section:

- Up to an eight character archive database name
- A path for the location of the archive databases
- A description for the archive databases

After completing the *Enter Archive Information* section, click on the <SAVE> button to save the archive databases or on the <CANCEL> button to cancel the archive databases.

A user may choose to delete an archive database if all of the records contained in it have been restored, it was created in error, or it was never utilized. To delete an archive database, select an archive database from the *Select an Archive Database* list and click on the <DELETE> button. Note that an archive database must be empty before it can be deleted.

To mark an archive database as the active archive database, select an archive database from the *Select an Archive Database* section. Click on the *Make this Database the Active Archive Database* check box and then click on the <SAVE> button to mark the selected archive database as the active archive database.


To exit the **Archive Setup** function, click on the <CANCEL> button.

Tracking Setup

The purpose of the **Tracking** function is to record significant events associated with the primary objects in the RAVEN-SB system: Resident; Facility; Employee; and Assessment. If the Tracking system is enabled,

then any events enabled for these objects will be captured and saved to the SWBTRACK.MDB database. This function enables the System Administrator to generate a history of events.

For example, as depicted in the screen below, each time an assessment is added, edited, or deleted, a record containing the Assessment ID, the Event ID, the Employee ID, and the time stamp of when the event occurred is created in the Assessment table of the SWBTRACK.MDB. In addition, each time an assessment is modified or inactivated, an appropriate event is generated.

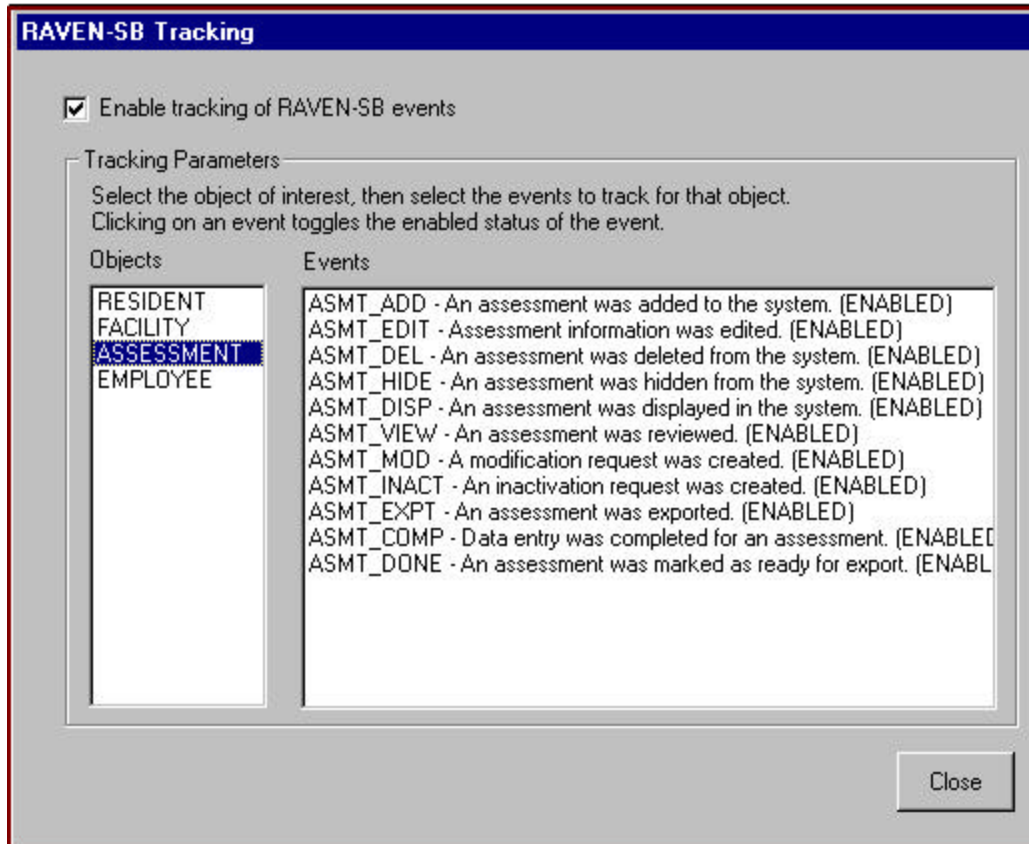


After a System Administrator selects the Tracking Setup option from the ADMINISTRATION menu item, the **RAVEN-SB Tracking** screen will be displayed. Click on the *Enable Tracking of RAVEN-SB Events* check box to enable the **Tracking** function. Select one of the objects, Resident, Facility, Assessment, or Employee, from the *Objects* list box to display a list of the possible events to track in the *Events* list box. Clicking on an event toggles the enabled status of that event. Click on an event to give it a status of "Enabled". Click on an enabled event to remove its enabled status. After defining the required Tracking events, click on the <CLOSE> button to exit the **Tracking** function.

NOTE:

To stop the tracking of RAVEN-SB events, yet maintain the tracking settings, click on the *Enable Tracking of RAVEN-SB Events* check box to disable the **Tracking** function.

Preferences

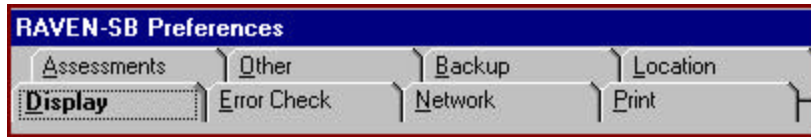

Display

Choose to Disable Color Tabs

When data entry is performed in RAVEN-SB, the main tabs are displayed in red and the subtabs are displayed in blue. RAVEN-SB gives the System Administrator the capability to turn off these colors and display grey tabs and subtabs. When grey tabs and subtabs are displayed, the name of the tab or subtab selected for data entry is bolded.

After selecting the Preferences option from the ADMINISTRATION menu and clicking on the Display tab, click on the *Disable Color Tabs* check box to disable the colored tabs/subtabs.

Error Check

Choose to Display Warning Messages When Error Checking Assessments

RAVEN-SB will display warning messages as part of the error checking process. These messages indicate whether assessments are in violation of the timing requirements as defined in the SB-MDS 1.0 Specifications. The messages do not prevent completion or submission of assessments. They are displayed for information purposes only.

The RAVEN-SB software defaults to display warning messages only when the <ERROR CHECK> button on the **Export** screen is clicked. After selecting the Preferences option from the ADMINISTRATION menu and clicking on the Error Check tab, click on the *Show Warning Messages when Error-Checking Assessments* check box to display warning messages during the **Data Entry** or **Import** processes.

NOTE:

Warning messages are always displayed when error checking from the **Export** screen.

Network

By default, a user has the option of accessing an assessment already marked as open for data entry. RAVEN-SB has the capability to allow only one user to access an assessment at a time. This option is useful for network environments.

After selecting the Preferences option from the ADMINISTRATION menu and clicking on the Network tab, click on the *Stop Users from Entering Open Assessments* check box. This selection will prevent multiple users from opening the same assessment.

Print

Specify a Path for the Adobe Acrobat Reader in Order to Print Assessment Forms

To print the RAVEN-SB assessment forms, the file/path name for the Adobe Acrobat Reader must be entered. The 32-bit Adobe Acrobat Reader executable name is ACRORD32.EXE. From the ADMINISTRATION menu item, select the Preferences option and then click on the Print tab and enter the file path/name into the *Adobe Acrobat Reader Executable File Path/Name* text box (e.g., C:\PROGRAMS FILES\ADOBE\ACROBAT5\READER\ACRORD32.EXE.).

NOTE:

A complete explanation of the **Print Assessment(s)** function is available in the System Setup section.

Enable Silent Printing

RAVEN-SB can prevent the printer dialog box from appearing when launching Adobe Acrobat Reader to print an assessment on the SB-MDS forms. After selecting the Preferences option from the ADMINISTRATION menu and clicking on the Print tab, check the *Enable Silent Printing* check box to enable this feature.

Assessments

Choose Radio Button or Check Box for Display of Option Pick One Variables

For variables that allow the user to select only one choice (e.g., Option Pick One), the default display uses radio buttons. For keyboard data entry that utilizes tabbing between variables, radio buttons automatically select the first item displayed in the list. Choosing check boxes instead of radio buttons for this variable type alters system behavior. This selection allows tabbing to and from Option Pick One-type variables without forcing the user to select an option. If an option is selected, RAVEN-SB will still enforce the rule that ensures only one option from the list of items can be selected.

After selecting the Preferences option from the ADMINISTRATION menu and clicking on the Assessments tab, click on the *Use Check Boxes for Variables Where You Select Only One Choice* radio button to change the display to check boxes.

Other

Specify Location Name and Telephone Number

By specifying a location name and telephone number for this installation of RAVEN-SB, employees who *log in* will be associated with this location. The User List menu option from the **RAVEN-SB Management** screen displays this information. After selecting the Preferences option from the ADMINISTRATION menu and clicking on the Other tab, type the location name and telephone number of the current installation in the *Location Name* and *Telephone* text boxes provided.

NOTE:

The location name and telephone number for a new employee will not be displayed until that employee logs into the system.

Enable Warning for Number of Assessments in Archive Database

By default, RAVEN-SB will not warn the user of the number of assessments contained in an archive database. However, RAVEN-SB gives the user the capability to display a warning when an archive database has reached a specific number of records.

After selecting the Preferences option from the ADMINISTRATION menu and clicking on the Other tab, click on the *Warn if an Archive Database Contains ____ or More Assessments* check box to display warning messages during the archiving process. The default value for the number of assessments on which to warn is 0. To change the default number, click in the box and enter the desired number, up to 99,999.

Backup

By specifying a path for the **Backup** function in RAVEN-SB and a Compression option, the collection databases in the RAVEN-SB system will be automatically backed up every time the system is exited normally. The path for the backup process cannot be the same as the RAVEN-SB installation path. For example, if RAVEN-SB is installed to C:\SWB, the backup path cannot be C:\SWB. If the path is left blank, the **Backup** function is disabled.

NOTE:

A complete explanation of the **Backup** function is available in the System Setup section.

Location

Specify Location of Databases

By default, RAVEN-SB looks for the data entry databases (SWB.MDB) in the application directory. To specify a different location for these files, select the Preferences option from the ADMINISTRATION menu item and click on the Location tab, then enter the desired path in the *Path* text box provided.

NOTE:

To use the new path, exit the system and restart the RAVEN-SB software.

Specify Default Paths for the RAVEN-SB Export, Import, and Print (from file) Functions

The user can optionally specify default locations for the following:

- Where export files are created
- Where import files are located
- Where print (from file) files are found

Note that for exports, the default output file name will be "TEMP.TXT". This file will be written to the export path, if specified. After selecting the Preferences option from the ADMINISTRATION menu item and clicking on the Location tab, enter the path locations for the desired features.

UNDER THE HELP MENU, THE USER CAN:

Access General Help and About Help

General Help provides RAVEN-SB system help regarding how to use and maintain the RAVEN-SB software. **About (RAVEN-SB)** identifies software and dictionary version information. From the **RAVEN-SB Management** screen or any of the **RAVEN-SB Data Entry** screens, click on the HELP menu and select either the General Help or the About options.


NOTE:

Context sensitive help is available from any window by pressing the <F1> key. Individual variable help is also available by right mouse clicking anywhere within a variable's frame.

SELECT /ADD A RESIDENT

The Resident Selection section enables the user to select the resident whose assessment history will be accessed. Residents are displayed in alphabetical order. Use the scroll bars on the *Resident* list box to view the available residents and click on the preferred row to select a resident. If there are more than 1,000 residents from which to choose, click on the <BROWSE NEXT 1000> button to display the next 1,000 residents.

The user may search the Resident database by selecting the Last Name, Social Security Number (SSN), Medicare Number, Medicaid Number, State Resident ID, or Facility Resident ID of a resident in the *Search for...* drop-down list box.

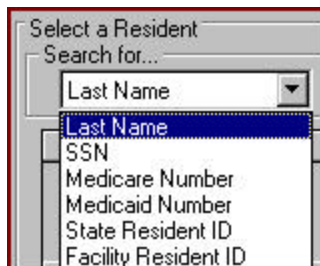

After selecting an item from the *Search for...* drop-down list box, click on the <SEARCH> button to begin the search. The closest match will be the first item displayed in the *Resident* list box. To remove the search, clear the = text box and click on the <SEARCH> button. The residents will then be displayed alphabetically for selection within the *Resident* list box.

The user may want to reduce the number of residents available from which to choose by filtering on Unit. This filtering is accomplished by entering the desired Unit in the *Filter for Unit...* text box. Click on the <FILTER> button to begin the filter. An alphabetical list of residents by unit will be displayed in the *Resident* list box. If the Unit list is extensive, the user can search this list by Last Name, Social Security Number (SSN), Medicare Number, Medicaid Number, State Resident ID, or Facility Resident ID. To remove the filter, clear the Unit field and click on the <FILTER> button. An alphabetical list of residents will then be displayed for selection within the *Resident* list box.

The user can also sort by a field (e.g., Last Name, First Name, MI, SSN, Unit, Medicare Number, Medicaid Number, State Resident ID, or Facility Resident ID) by clicking on the column header. Clicking once on the

column header will sort the list in ascending order. Clicking on the column header again will sort the list in descending order.

The <ADD RESIDENT> button allows the Clinical Supervisor or System Administrator to add a new resident record. The user must complete the Last Name, Gender, Birthdate, and Race/Ethnicity before the record can be saved. The remaining fields are optional.

To edit, hide, or delete an existing resident record, the Maintain Resident Database option from the ADMINISTRATION menu item must be selected. Refer to the Maintain Resident Database section for a detailed explanation of how to add/edit/hide/delete a resident.

ARCHIVE A RESIDENT

The purpose of the **Archive** feature is to help maintain the size of the system database SWB.MDB. The **Archive** function allows the System Administrator to archive a resident and his or her associated assessments when the last completed assessment was a discharge assessment and all of the assessments for that resident have a status of "Exported". For instance, if a resident has been discharged due to death, discharged with no intention of returning, or a specified amount of time has passed since his/her last discharge assessment, he/she and all of his/her assessments can be archived for storage purposes, thus reducing the size of the system database SWB.MDB.

To begin the **Archive** function, select a resident whose last completed assessment was a discharge assessment and whose assessments have all been exported and click on the <ARCHIVE RESIDENT> button. RAVEN-SB will archive that resident's assessments to the active archive database. The System Administrator also has the capability to create an archive database and specify the location of that archive database. Click on the ADMINISTRATION menu item and select the Archive Setup option to display the **Archive Information** screen. See the Archive Setup section under the Data Entry/Select Case Functions section for a complete explanation of creating an archive database.

NOTE:

To view archived resident records on the **Resident Information** screen, click on the ADMINISTRATION menu item, select the Maintain Resident Database option, and click on the *View Archived Residents Also* check box. To view archived resident records from the **RAVEN-SB Management** screen, click on the FILE menu item and select the Open Archive option.

It is recommended that the system database, SWB.MDB, be repaired and compacted after several archives have been performed. In order to complete the **Repair/Compact** process, all users must exit the RAVEN-SB software. See Appendix A for a complete description of executing the **Repair/Compact** function.

RESTORE A PATIENT

The **Restore** function allows the System Administrator to reinstate a resident who has been archived and removed from the system database SWB.MDB. To begin the **Restore** function, select the Open Archive option from the FILE menu to display the residents who have been archived in the default archive database or in the archive database specified on the **Archive Information** screen. All archived residents will be displayed on the **RAVEN-SB Management** screen. Select a resident and click on the <RESTORE RESIDENT> button to reinstate that resident. To close the archive database and display the active residents in the system database, SWB.MDB, reselect the Open Archive option from the FILE menu.

While in the **Restore** mode of RAVEN-SB, the user has limited functionality. These features include: switching between the **Archive** and **Restore** functions; the **Reports** function; the **Log Off/On** function; accessing **Help**; viewing and printing assessments; and RUGs.

SELECT /ADD AN ASSESSMENT

The Assessment Selection section allows the user to view the history of previous assessments and their status for the currently selected resident. For each assessment, the date, assessment type, status, last user to access the assessment, and open status (an asterisk [*] indicates an assessment is currently open for data entry) are displayed.

NOTE:

An assessment is "Open" once it is accessed for data entry. The open status only occurs when an assessment is accessed in a network installation of RAVEN-SB. The *Open* box will be empty when a user is working in a stand-alone configuration.

An assessment is selected by scrolling through the list and clicking on the preferred assessment. The user can access the selected assessment to perform data entry by clicking on the <ENTER DATA> button or by double-clicking on the assessment. The user may run RUGs (both discussed in this manual) or print the highlighted assessment by clicking on the appropriate button.

The Add Assessment section also enables the user to initiate a new assessment for the selected resident. Click on the <ADD ASSESSMENT> button to access the **Reasons for Assessment -- Add Assessment** screen (described below) before the assessment is added to the resident.

Note the following:

- Newly created assessments have a status of "New"
- Edited assessments have a status of "In Use"
- Exported assessments have a status of "Exported"
- Closed assessments (data entry completed) have a status of "Export Ready"
- Hidden assessments have a status of "Hidden"

Hidden assessments are only visible when "View Hidden Assessments Also" is checked. This feature is automatically enabled when an assessment is hidden.

Reason for Assessment

When the user clicks on the <ADD ASSESSMENT> button, the **Reasons for Assessment -- Add Assessment** screen is displayed.

Reasons For Assessment - Add Assessment

Enter reasons for assessment as taken from item 11:

11. REASONS FOR ASSESSMENT

<p>a. Primary Reasons for Assessment</p> <p><input checked="" type="radio"/> 00. PPS assessment for Medicare Payment</p> <p><input type="radio"/> 06. Discharged-Return Not Anticipated</p> <p><input type="radio"/> 07. Discharged-Return Anticipated</p> <p><input type="radio"/> 09. Reentry</p> <p><input type="radio"/> 11. Assessment-Not for Medicare payment</p>	<p>b. PPS Scheduled Assessments</p> <p><input type="radio"/> 1. 5-day <input type="radio"/> 7. 14-day</p> <p><input type="radio"/> 2. 30-day <input type="radio"/> 9. Other</p> <p><input type="radio"/> 3. 60-day</p> <p><input type="radio"/> 4. 90-day</p> <p><input type="radio"/> 5. Readmission/Return</p> <p>Clear Value</p>
--	---

☐ c. OMRA Assessment
☐ d. Clinical Change Assessment
☐ e. State-Required Assessment
☐ f. Assessment Needed for Other Reasons (e.g., HMOs, MSP, sanction situations, etc.)

Add Cancel

The selection made in Section 11a (Primary Reason for Assessment) will determine the options available in each subsequent section. After the selections have been made, click on the <ADD> button to save those choices and return to the **RAVEN-SB Management** screen or click on the <CANCEL> button to cancel the assessment and return to the **RAVEN-SB Management** screen.

NOTE:

If the RAVEN-SB Sequencing of Assessments option is selected, the user will be restricted to adding assessments that can sequentially follow the assessment highlighted on the **RAVEN-SB Management** screen.

Correct an Assessment

Correction Type

Specify the type of correction to be made:

☒ Modify the assessment
☐ Inactivate the assessment

OK Cancel

Once an assessment has been Exported only a System Administrator has access to it. Listed below are the two types of corrections that can be performed:

- Assessment was Submitted to the National Database and was Accepted. Correction to Fields is Necessary.** By selecting this option, a copy of the assessment record is created with an incremented correction number, and the user is taken directly into data entry for that assessment. After the necessary corrections are made, the assessment can be completed, exported, and re-submitted to the State System.
- Assessment was Submitted to the National Database and was Accepted. Inactivation of the Assessment is Necessary.** By selecting this option, a copy of the assessment record is created. The user is taken directly into data entry for that assessment, however the entire assessment is view-only.

After viewing the assessment, the user can exit the assessment. The correction number on the **RAVEN-SB Management** screen will be 99.

NOTE:

After exporting an inactivation record, a new assessment containing the data from the inactivated assessment can be created. Select the inactivated assessment from the **RAVEN-SB Management** screen and click on the <CORRECT ASSESSMENT> button. A confirmation box will appear. Click on the <OK> button to create the assessment prepopulated with the data from the inactivated assessment. The correction number for this assessment will be 00.

SPECIAL FEATURES (RAVEN-SB MANAGEMENT SCREEN)

View Only (Data Entry Not Allowed)

The View Only (Data Entry Not Allowed) option is only visible to users with Clinical Supervisor or System Administrator rights. If checked, this option allows the Clinical Supervisor or System Administrator to view an assessment in read-only mode. This option is useful in reviewing the work of other users while ensuring that no data is accidentally edited.

View Hidden Assessments Also

The View Hidden Assessments Also option is available only to users with Clinical Supervisor or System Administrator rights. If checked, this option enables the user to see and access previously hidden assessments (i.e., assessments saved in the database but not available for viewing by users with Data Entry privileges only).

EXIT

The user may exit the RAVEN-SB software by clicking on the <EXIT> button or by selecting the Exit option from the FILE menu item.

RAVEN-SB DATA ENTRY WINDOW

The **Data Entry** screen consists of the following components:

- **Title Bar.** Displays the RAVEN-SB function **Data Entry**.
- **Menu Bar.** Displays the list of functions that can be performed on the **Data Entry** window. Each of these functions contains additional options that can be displayed by using the mouse.
- **Data Entry Button Bar.** Provides the available **Data Entry** functions. Currently, the following functions are available: **Enter Notes** (for an assessment), and **Pause Data Entry** (pause abstraction time for data entry).
- **Data Entry Top.** Displays the Last Name, First Name, Social Security Number (SSN), and Birth Date of the selected Resident at the top of the **Data Entry** window.
- **Tabs/Subtabs.** Displays the names of the screens for the selected set of screens. Each tab can have up to five subtabs associated with it. When selecting a tab or subtab, the **Data Entry** window of the selected screen or subscreen will open in the space below it.
- **Data Entry Window.** Contains variables that belong to the selected screen or subscreen. On this screen, data entry can be performed for each variable.


ENTER DATA FOR A VARIABLE

Perform data entry after selecting a Resident record and an Assessment form. See previous sections, Selecting/Adding a Resident and Selecting/Adding an Assessment for further information. Data entry can be exited at any time by clicking on the FILE menu item and selecting the Exit and Validate Assessment option.

Using the mouse is the most efficient way to move around the data entry system. Tabbing is also available. Since Microsoft™ has restrictions pertaining to variables with radio buttons, the mouse must be used to move between variables of this type (**Option (Pick One)**).

To begin data entry, point to the preferred tab or subtab and click on it to open the corresponding screen. Follow the descriptions below to enter data for the various types of variables.

Option (Pick One)

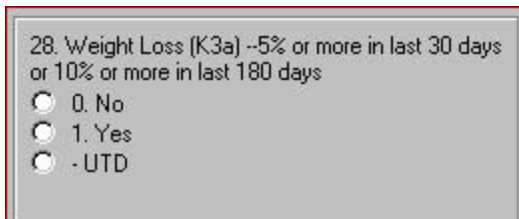

28. Weight Loss (K3a) --5% or more in last 30 days
or 10% or more in last 180 days


☐ 0. No

☐ 1. Yes

☐ - UTD

The variable type **Option (Pick One)** includes the options from which only one option can be chosen. Choose an item by selecting the radio button next to the item. Clicking on the <RIGHT MOUSE> button and selecting the Clear Variable option, clears variables with radio buttons.

Option (Pick One or More)

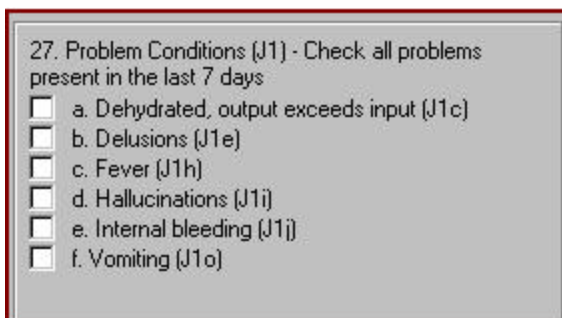

27. Problem Conditions (J1) - Check all problems
present in the last 7 days

☐ a. Dehydrated, output exceeds input (J1c)

☐ b. Delusions (J1e)

☐ c. Fever (J1h)

☐ d. Hallucinations (J1i)

☐ e. Internal bleeding (J1j)


☐ f. Vomiting (J1o)

The variable type **Option (Pick One or More)** includes options from which one or more items can be chosen. Choose the item(s) by marking the check box(es) next to the desired item.

NOTE:

Check boxes have three possible states: checked, unchecked, and checked/grayed. Checked indicates that the option is true. Unchecked indicates that the option is false. Checked/grayed indicates that the option is Unable to Determine (UTD). To set the option to UTD, click the <RIGHT MOUSE> button on the check box and select the Set as UTD option.

Date

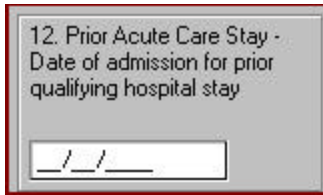


12. Prior Acute Care Stay -
Date of admission for prior
qualifying hospital stay

__/__/____

The variable type **Date** accepts data entered in the date format that has been specified in the design (the default format for date is MM/DD/YYYY).

If the date cannot be determined, press the hot key combination **<SHIFT>** and **<?>** (shift and question mark keys). A dash will be displayed.

Number

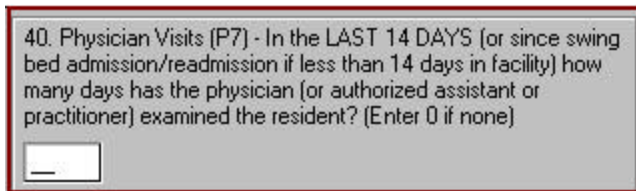

40. Physician Visits (P7) - In the LAST 14 DAYS (or since swing
bed admission/readmission if less than 14 days in facility) how
many days has the physician (or authorized assistant or
practitioner) examined the resident? (Enter 0 if none)

__

The variable type **Number** accepts data entered in numeric format, including one decimal point, if applicable.

If the value cannot be determined, press the hot key combination **<SHIFT>** and **<?>**. A dash will be displayed.

String



c. (Last)

GERSHWIN

The variable type **String** accepts up to 99 characters in any alphanumeric combination.

If the value cannot be determined, press the hot key combination **<SHIFT>** and **<?>**. A dash will be displayed.

SPECIAL FEATURES (RAVEN-SB DATA ENTRY SCREEN)

Note Box

To enter a general note about a variable during data entry, select the <NOTES> button. The box allows the user to enter a note containing up to 32,000 alphanumeric characters. Since all the notes for a given assessment are stored in one file, copying the variable screen title before writing a note provides a good reference point to the variable. In order to paste the screen title of the current variable at the bottom of the note, select the <COPY> button. Click on the <CLOSE> button to save the notes and return to the **Data Entry** window.

Pause Button

When the user selects the <PAUSE ABSTRACTION> button, the internal RAVEN-SB clock that keeps track of the amount of time spent on an assessment is stopped.

NOTE:

When in **View Only** mode, abstraction time does not accumulate.

Clear Value on a Variable Radio Button

To clear variables with radio buttons, click on the <RIGHT MOUSE> button and select the Clear Variable option.

View Help

To display the Help available for each variable, click on the <RIGHT MOUSE> button and select the View Help option. For variables that do not have any other options associated with them (e.g., the UTD option), the help will be displayed immediately.

Print Screen

Selecting the Print Screen option from the FILE menu item allows the user to print the current screen. If the screen does not print properly (possibly due to the default printer settings), another process, screen capture, is available.

Several steps need to be followed when capturing a screen for printing. To copy a screen, press on the <PRINT SCRIN> button or press on the <ALT> + <PRINT SCRIN> buttons to capture only the active window. After opening a word processing package, such as Word, or a graphics package, such as Paint, paste the screen into the document by pressing on the <CTRL> + <V> buttons. The screen also can be copied into a document by using the paste function available in the application being used.

Edit Facility Data

During data entry, the user will have access to the **Facility Information** screen by selecting the Edit Facility Data option from the ADMINISTRATION menu item. This feature allows the user to correct the Facility Information for an assessment.

For example, if incorrect information is added during the setup process for the facility, the user can edit this information during the data entry process. This edit will correct the information for this assessment and future assessments. Prior assessments are not automatically corrected. It is necessary to revisit those assessments to update their facility information.

After returning to the **RAVEN-SB Management** screen, select each assessment (by resident) that contains the incorrect information and click on the <ENTER DATA> button. A message indicating that the data have been changed for the facility will appear and the user has the option to update the current assessment.

NOTE:

Only the Facility and Agent Information can be edited on the **Facility Information** screen during **Data Entry** mode. The Facility/State Options cannot be edited while in **Data Entry** mode.

Edit Resident Data

During data entry, the user will have access to the **Resident Information** screen by selecting the Edit Patient Data option from the ADMINISTRATION menu item. This feature allows the user to correct the Resident Information for an assessment.

For example, if incorrect information is added during the setup process for a resident, the user can edit this information during the data entry process. This edit will correct the information for this assessment and future assessments.

Prior assessments are not automatically corrected. It is necessary to revisit those assessments to update their resident information. After returning to the **RAVEN-SB Management** screen, select each assessment that contains the incorrect information for the resident and click on the <ENTER DATA> button. A message indicating that the data have been changed for the resident will appear and the user has the option to update the current assessment.

NOTE:

When editing the **Resident Information** screen during data entry, only the Resident Information portion of the screen will be visible. The user cannot search for or select another resident.

CALCULATE/VIEW RUGs (IF APPLICABLE)

Calculate RUGs

RUG Values:

43MDCR. Medicare case mix group Version

Calculation Method

CMI Set

43STATE. State case mix group Version

Calculation Method

CMI Set

NOTE: RUGs are stored with this assessment.

RUG calculations are automatically performed for an assessment when the data entry for the assessment has been completed and its status is made "Export Ready". RUGs can be calculated for all applicable assessments. The RUGs can be recalculated at any time during data entry. Click on the RUGs menu item in data entry to recalculate the RUGs. The **Calculate RUGs** screen will be displayed. Note, however, that the RUG values are not saved until the assessment is locked.

RUGs can also be viewed for an assessment from the **RAVEN-SB Management** screen using the <RUGS> button. If the assessment is "Export Ready" or "Exported", the RUGs have already been calculated and will merely be displayed in the **View RUGs** screen. RUGs will be recalculated and displayed for applicable assessments with other statuses. Note that the calculation method and CMI set are not displayed, since these parameters may have been modified by the user since the RUGs were calculated and stored for the assessment.

See the Maintain Facility Database section in the System Setup chapter for a complete description of how to setup the RUG calculations.

SB_SUB_REQ

SB_SUB_REQ

Select the appropriate SB_SUB_REQ value for this assessment:

☒ 1 - If the patient IS in a Medicare- or Medicaid-certified swing bed. The Swing Bed-MDS event (assessment reference date, discharge date, or reentry date) DID occur while the patient was in a Medicare- or Medicaid-certified swing bed.

☐ 0 - If the patient IS NOT in a Medicare- or Medicaid-certified swing bed. The Swing Bed-MDS event (assessment reference date, discharge date, or reentry date) occurred while the patient was in a bed other than a certified swing bed.


Save Cancel

Click on the SB_SUB_REQ menu item to display the **SB_SUB_REQ** screen. All assessments are required to have a SB_SUB_REQ value.

ASSESSMENT VALIDATION AFTER EXIT

When the user attempts to exit data entry, the system will check all of the rules in place. As the rules are being processed, the status of the validation will be displayed in the center of the screen. One such rule is the mandatory variable rule. All mandatory variables must have valid data values. If a mandatory variable is left blank (i.e., data are not entered), the user will be warned that this variable is mandatory (e.g., ERROR: A value is required for this field.) The user will not be able to mark an assessment as complete until all mandatory variables have been entered.

In addition, any assessment validation warnings will be displayed. Since these warnings are for informational purposes and are not considered errors, they are not required to be resolved before an assessment can be marked as complete.

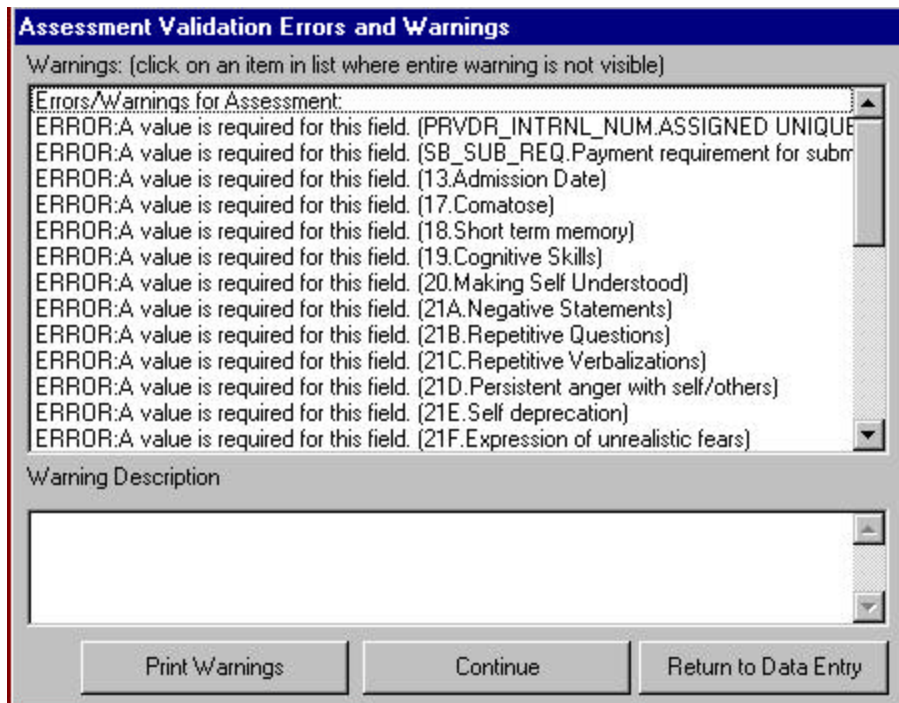


Three options are available on the **Assessment Validation Errors and Warnings** screen. Clicking on the <PRINT WARNINGS> button will print the first 30 errors/warnings listed on the screen. Selecting the <CONTINUE> button will take the user to the **Assessment Status** screen where the assessment can be exited by marking it as complete or in use. Choosing the <RETURN TO DATA ENTRY> button will allow the user to continue the data entry process.

NOTE:

To enable the display of warning messages for timing requirements, the System Administrator must select the Error Check tab from the PREFERENCES menu item and then click on the *Show Warning Messages when Error-Checking Assessments* check box.

ASSESSMENT STATUS WINDOW

On the **Assessment Status** screen, the status of an assessment must be marked before the assessment can be exited.

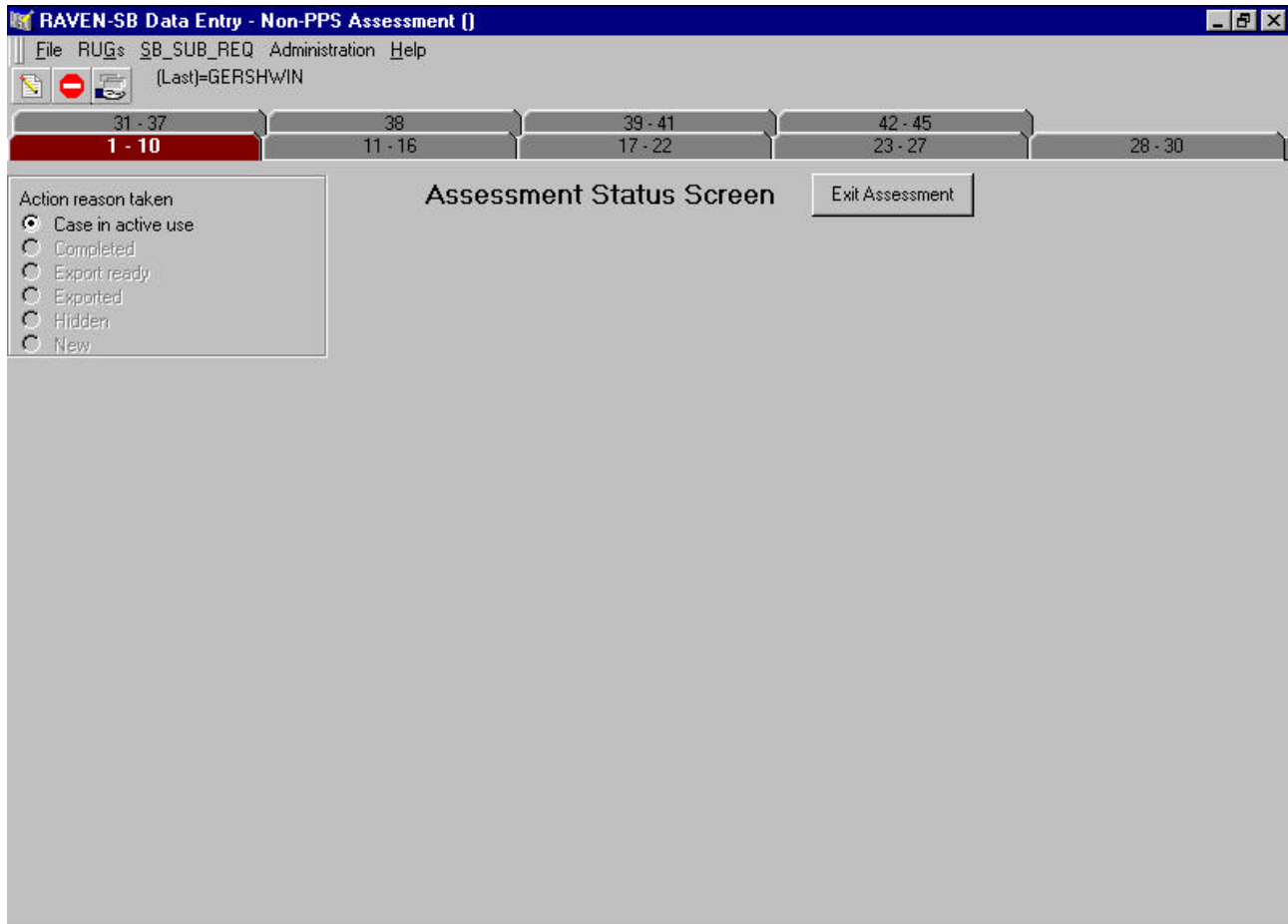

The various options that can be selected when exiting the **Assessment Status** screen include:

- Assessment in active use -- An assessment has been edited
- Completed -- Data entry has been completed for an assessment

The options that only can be viewed when exiting the **Assessment Status** screen include:

- New -- A newly created assessment
- Exported -- An exported assessment
- Hidden -- A hidden assessment
- Export Ready -- An assessment with a status of closed (data entry completed)

At a minimum, the user will be required to mark the assessment as "Assessment in active use" or "Completed". These two options are the default options on the **Assessment Status** screen.

In addition, the other options available on the **Assessment Status** screen, New, Exported, Hidden, and Export Ready, will automatically be selected depending on the type and status of the Assessment.

To mark the status, select the best choice and press the <EXIT ASSESSMENT> button.

NOTE:

Once an assessment has been marked "Export Ready" or "Exported", only users with Clinical Supervisor or System Administrator rights are permitted to reopen it.

RAVEN-SB DATA EXPORT

EXPORTING ASSESSMENTS

To initiate the **Export** function, select either the Export option from the OPTIONS menu item or the <EXPORT> button on the **RAVEN-SB Management** screen.

RAVEN-SB Export

Export Details

☒ New Export (Assessments Marked as Export Ready)
☐ Previous Export (Assessments Marked as Export Ready or Exported)

☒ Production Submission
☐ Test Submission

Export Path and File Name (including extension): C:\SWB\

Export Description (optional):

Previous Exports:

- 34 (11/01/2001, C:\JUNK.TXT, 2 assessments)
- 33 (11/01/2001, C:\JUNK.TXT, 2 assessments)
- 32 (11/01/2001, C:\JUNK.TXT, 2 assessments)

Select Assessments to be Included in Export

#	Name	11A	11B	Date	43MDCR	43STATE	Status
---	------	-----	-----	------	--------	---------	--------

☐ Select All ☐ Clear All

Print Assessment List

Error Check Export Close

The RAVEN-SB **Export** function enables data for "Export Ready" assessments to be output from the RAVEN-SB database into the Data Record Layout format (i.e., 1,814-byte string) suitable for submission to the State. Note that each data field in each exported assessment has been error-checked in accordance with the rules specified in the RAVEN-SB Data Record Layout.

Following is a description of the **Export** function.

- Select the type of submission -- New Export or Previous Export. Previous Exports are described in the next section.
- After selecting the New Export radio button, a description of the file may be entered into the *Export Description* text box. If no description is entered, a default description will be generated stating the number of assessments exported. The assessments will be displayed in the *Select Assessment to be Included in Export* text box. Select the assessment(s) to be included with this export from the *Select Assessments to be Included in Export* text box. To select one assessment, click on the desired assessment. To select more than one assessment, use the <SHIFT> or <CTRL> keys. To include all of the displayed assessments for export, click on the Select All radio button. Click on the <CLEAR ALL> button to deselect all highlighted assessment.
- Enter the path/file name into the *Export File Path/Name* text box. Either type the information in the text box or click on the <FILE FOLDER> icon, select a path and enter a file name. The data will be saved in text file format. The user may want the file name to correspond to the description used in the *Export Description* text box.

EXAMPLE:

For easy reference, the user may wish to create a directory called RAVEN-SB with subdirectories for the year/month. The saved file names could reference the weeks of the month (i.e., C:\SWB\2001\OCTOBER\WEEK1.TXT). This structure will help to organize the export files for submission to the State System.

- Select the type of submission -- Production or Test. An assessment receives the status "Exported" once it has been exported using the Production Submission option. Note that the user has the capability to generate a test submission (i.e., create a test export file) without changing the status of the assessment records. This structure allows the agency to validate the submission format with the State System.
- Click on the <ERROR CHECK> button to check the highlighted assessments for errors. This action should result in warning messages for timing violations only since each assessment is "Export Ready". The user has the option to print any warning messages that appear.
- Click on the <EXPORT> button to complete the export process. A message will be displayed stating that the export is complete and the file is ready to submit to the State. Click on the <OK> button to close the message box. After the export has been completed, the current date, file and path name, and a description of the exported file will be displayed in the *Previous Exports* text box. Click on the <CLOSE> button to exit the **Export** function and return to the **RAVEN-SB Management** window.

NOTE:

Some assessments may not be included in the submission file, however the status of those assessments will change from Export Ready to Exported. The excluded assessments include: (1) Assessments that have SB_SUB_REQ=0.


RE-EXPORTING ASSESSMENTS

It also is possible to resubmit previously exported assessments, should an error be detected by the State System. To create a resubmission file, the user selects the Previous Export radio button and is presented with a list of previous submissions in the *Previous Exports* list box. After selecting a previous submission, the user may then generate a new export file that contains data for each of the assessments displayed in the *Select Assessments to be Included in Export* text box.

RAVEN-SB DATA IMPORT

SELECTING IMPORT OPTIONS

The System Administrator has the capability to import files. Import files must utilize the same format as export files (i.e., they must contain a header record, assessments in the 1,814-byte string format specified in the SB-MDS Data Record Layout, and a trailer record). On the **RAVEN-SB Management** screen, the user may click on the <IMPORT> button or select the **Import** function from the OPTIONS menu item to display the **RAVEN-SB Import Assessments** screen.

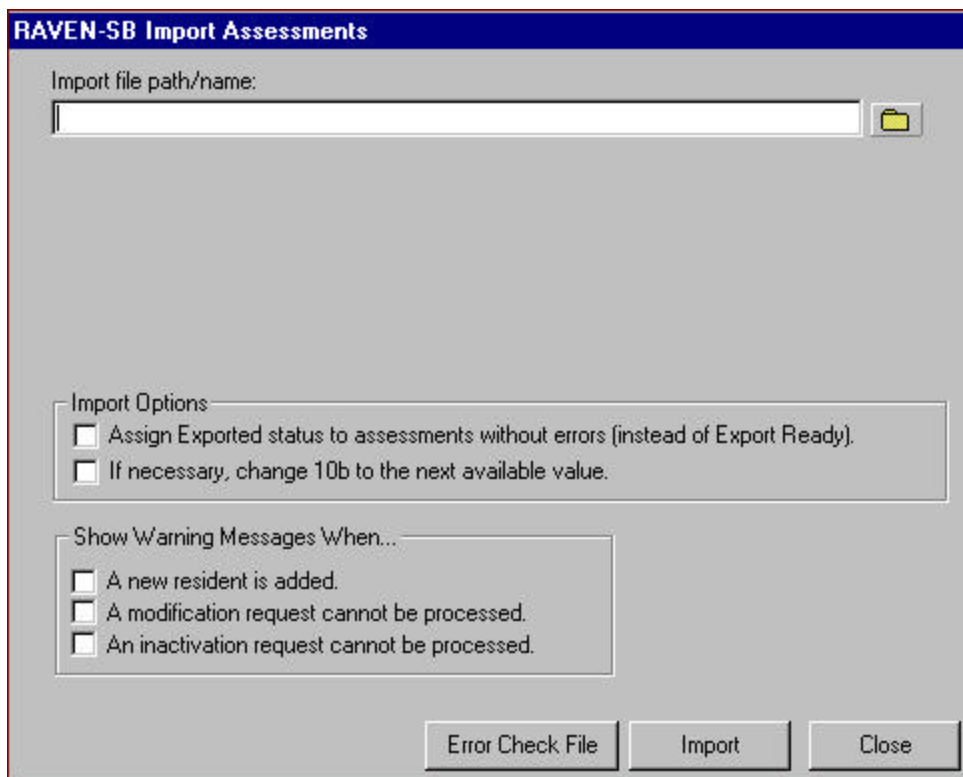

The screenshot shows a dialog box titled "RAVEN-SB Import Assessments". It features a text input field for "Import file path/name:" with a folder icon to its right. Below this, there are two sections of options. The first section, "Import Options", contains two checkboxes: "Assign Exported status to assessments without errors (instead of Export Ready)." and "If necessary, change 10b to the next available value." The second section, "Show Warning Messages When...", contains three checkboxes: "A new resident is added.", "A modification request cannot be processed.", and "An inactivation request cannot be processed." At the bottom of the dialog are three buttons: "Error Check File", "Import", and "Close".

The user must enter an Import file path/name for the file to be imported or use the default file path/name. The default file path is the directory where the RAVEN-SB software was initially installed with a file name of TEMP.TXT.

Two Import options are available when importing files. They include:

- **Assign Exported Status to Assessments without Errors (Instead of Export Ready).** This feature allows users to specify how completed assessments with no errors are handled by the **Import** function. By default, these assessments are given the status of "Export Ready". However, if this feature is checked, the assessment will be given the status of "Exported".

- **If Necessary, Change 10B to the Next Available Value.** If a modification or inactivation request is valid except for the 10B value, RAVEN-SB can automatically assign the next available value to the 10B variable and proceed with the import of the modification or inactivation.

SELECTING TO SHOW WARNING MESSAGES

During the Import process, the user has the capability to display several types of warning messages. These include: a new resident is added; a modification request cannot be processed; and an inactivation request cannot be processed. Select the desired warnings to be displayed during import by clicking in one or more of the warning check boxes.

For uninterrupted processing, leave these options unchecked.

PERFORMING ERROR CHECKS ON AN IMPORT FILE

Each record in the import file is processed in accordance with the rules specified in the SB-MDS Data Record Layout. The user may check the file for errors by clicking on the <ERROR CHECK FILE> button.

If errors are encountered, a list box describing the line number, field name, and type of each error encountered, as well as the number of records processed, will appear. The user may then click on the <PRINT> button to print the error list or the user may click on the <CLOSE> button to return to the **RAVEN-SB Import Assessments** screen. If no errors are detected, a message box appears to announce that the error check is complete and to display the number of records processed.

NOTE:

This function does not import data from any other type of file format.

IMPORTING RAVEN-SB DATA

The user may proceed with the import by clicking on the <IMPORT> button. This action will process all records in the designated import file. The following paragraphs describe how the **Import** function handles facility, resident, and assessment data.

- **Facility Data**

The Facility information for the header record in the import file is compared with the existing facility information in the SWB database. If the information does not match, the user is given the option to update the facility database.

- **Resident Data**

The Resident information for each assessment in the import file is matched against the data in the existing resident database. The software compares the following fields: First Name, Last Name, Social Security Number, and Gender. If a match is not found, a new entry will be added to the resident database.

By default, a message will not be displayed when a new resident is added; however, the user can override this feature by checking the box labeled "A new resident is added". This feature is useful when importing

assessments for a small number of residents who do not have previous resident database information in RAVEN-SB.

- **Assessment Data**

Imported assessments may either have a status of "In Use", "Export Ready" or "Exported". In order for an imported assessment to be "Export Ready" or "Exported", depending on the selected option, the assessment must pass all error checking. All other assessments will be assigned the "In Use" status.

When a modification assessment is imported, the RAVEN-SB software will verify whether or not the prior version of the assessment exists. If the prior assessment is not found, RAVEN-SB will always import the modification assessment. However, if the prior assessment is found, RAVEN-SB will only import the modification assessment if the prior assessment has a status of "Exported".

When importing an inactivation assessment, the RAVEN-SB software must locate the prior version of the assessment and the status of the prior assessment must be "Exported". If RAVEN-SB cannot locate the prior version of the assessment or the status of the prior assessment is not "Exported", the inactivation assessment will not be imported. If RAVEN-SB can locate the prior version of the assessment and the status of the prior assessment is "Exported", the import of the modification assessment will continue.

RAVEN-SB HELP

GENERAL/ABOUT HELP

The **General Help** provides information about how to use and maintain the RAVEN-SB software. The **About (RAVEN-SB)** identifies software and dictionary version information. From the **RAVEN-SB Management** screen, or any of the **RAVEN-SB Data Entry** screens, click on the HELP menu and select either the General Help or About options.

ADDITIONAL HELP

Context sensitive help is available from any window by pressing the <F1> key. Individual variable help is also available by right mouse clicking anywhere within a variable's frame.

REPORTS

SELECT A REPORT TO VIEW/PRINT

RAVEN-SB provides the user with a report that can be used to review the quality of the abstracted data as well as a report that tracks various events within the RAVEN-SB software.

NOTE:

The Event Tracking Report is only available to a System Administrator.

When the **Reports** function is selected, the user will be presented with a **Report Selection** window where the type of report can be specified. After a report type is selected, the report criteria for that report will be displayed.

For each of these reports, a temporary file is always produced. Therefore, it is necessary to specify the location and the name of the file in order to change the default file location and the file name provided by RAVEN-SB.

ASSESSMENT DATA ENTRY REPORT

This report summarizes the data entered for an assessment(s) for quality control/review purposes.

View/Print Assessment Data Entry Report

Select a report: Assessment Data Entry Report
Event Tracking Report

Search For Assessments By:
☒ Resident ☐ Status

Select A Status:
 Search for... In Use Search

Select A Resident
 Search for... Last Name = Search

Filter for...
 Unit = Filter

Last Name	First Name	MI	SSN	Unit
GERSHWIN	GEORGE		894-35-4986	5

Select Assessment(s) Below To Report

OK Cancel

After selecting the Assessment Data Entry Report, select either the Resident option or the Status Option to determine the assessments to be used in the report.

REPORT ON AN ASSESSMENT(S) USING THE RESIDENT FUNCTION

- **Search for an Assessment by Resident**

Click on the <RESIDENT> radio button and from the *Search for...* drop-down list box and select the desired criteria on which to search. Follow the instructions discussed in the RAVEN-SB Data Entry section under Select/Add a Resident. Click on the <SEARCH> button to perform the search and display the results in the *Select a Resident* list box. Resident search options include:

- **Last Name.** A list of residents will be displayed by Last Name in the *Select a Resident* list box.
- **SSN.** A list of residents will be displayed by Social Security Number in the *Select a Resident* list box.
- **Medicare Number.** A list of residents will be displayed by Medicare Number in the *Select a Resident* list box.
- **Medicaid Number.** A list of residents will be displayed by Medicaid Number in the *Select a Resident* list box.
- **State Resident ID.** A list of residents will be displayed by State Resident ID in the *Select a Resident* list box.

- **Facility Resident ID.** A list of residents will be displayed by Facility Resident ID in the *Select a Resident* list box.

NOTE:

To select assessments that are new or deleted, it is necessary to search by Resident.

After selecting a resident, the associated assessments will be displayed in the *Select Assessment(s) Below to Report* list box on which to report. Selecting one or more assessments on which to report, click on the <OK> button to produce the report(s) or the <CANCEL> button to return to the **View/Print Assessment Data Entry Report** screen. To exit the **Reports** function, click on the <CANCEL> button.

NOTE:

A report(s) can take some time to produce depending on the size of the data entry system. Also, if a large number of assessments exist in the system, it is not advisable to select all of the assessments at once because of the amount of time it takes to produce a report(s).

REPORT ON AN ASSESSMENT(S) USING THE STATUS FUNCTION

- **Search for an Assessment by Status**

Click on the <STATUS> radio button and from the drop-down list select the status on which to report. Click on the <SEARCH> button to perform the search and display the results in the *Select Assessment(s) Below to Report* list box. Status search options include:

- **In Use.** An alphabetical list of residents with assessments coded "In Use" will be displayed.
- **Export Ready.** An alphabetical list of residents with assessments coded "Export Ready" will be displayed.
- **Exported.** An alphabetical list of residents with assessments coded "Exported" will be displayed.
- **Hidden.** An alphabetical list of residents with assessments coded "Hidden" will be displayed.

After selecting one or more assessments on which to report, click on the <OK> button to produce the report(s) or the <CANCEL> button to return to the **View/Print Assessment Data Entry Report** screen. To exit the **Reports** function, click on the <CANCEL> button.

NOTE:

A report(s) can take some time to produce depending on the size of the data entry system. Also, if a large number of assessments exist in the system, it is not advisable to select all of the assessments at once because of the amount of time it takes to produce a report(s).

EVENT TRACKING REPORT

The purpose of the **Tracking** function is to record significant events associated with the primary objects in the RAVEN-SB System: Resident; Facility; Assessment; and Employee. If the Tracking function is enabled, then any events enabled for the primary objects will be captured and saved to the SWBTRACK.MDB database. The Event Tracking Report summarizes the selected events for each object and enables the System Administrator to generate a history of these events.

View/Print Assessment Data Entry Report

Select a report: Assessment Data Entry Report
Event Tracking Report

Select an object: RESIDENT Find Object Instances

The following object instances have events in the tracking database.
Select an instance from the list, then click OK to view its event history.

FIELDS: Last Name; BirthDate; SSN; Unit

GERSHWIN;05/05/1950;894354986;5;

OK Cancel

After selecting the Event Tracking Report, click on the drop-down menu list to select an object on which to report, and click on the <OK> button.

Select an object: RESIDENT

The following object instances have events in the tracking database.
Select an instance from the list, then click OK to view its event history.

RESIDENT
FACILITY
ASSESSMENT
EMPLOYEE

Once an object has been selected, click on the <FIND OBJECT INSTANCES> button to display the instances associated with that object. Select an instance, then click on the <OK> button to display its event history.

Following is a list of the events available for each of the objects in the RAVEN-SB system.

- **Resident**

RES_ADD -- A resident was added to the system.

RES_DEL -- A resident was deleted from the system.

RES_EDIT -- Resident information was edited.

RES_ARCH -- A resident was archived.

RES_HIDE -- A resident was hidden from the system.

RES_DISP -- A resident was displayed in the system.

RES_REST -- A resident was restored from the archive.

- **Facility**

FAC_EDIT -- Facility information was edited.

AGT_EDIT -- Agent information was edited.

RUG_EDIT -- RUG setup information was edited.

CMI_EDIT -- CMI setup information was edited.

CMI_ADD -- A CMI set was added to the system.

CMI_DEL -- A CMI set was deleted from the system.

- **Assessment**

ASMT_ADD -- An assessment was added to the system.

ASMT_EDIT -- Assessment information was edited.

ASMT_DEL -- An assessment was deleted from the system.

ASMT_HIDE -- An assessment was hidden from the system.

ASMT_DISP -- An assessment was displayed in the system.

ASMT_VIEW -- An assessment was reviewed.

ASMT_MOD -- A modification request was created.

ASMT_INACT -- An inactivation request was created.

ASMT_EXPT -- An assessment was exported.

ASMT_COMP -- Data entry was completed for an assessment.

ASMT_DONE -- An assessment was marked as ready for export.

- **Employee**

EMP_ADD -- An employee was added to the system.

EMP_DEL -- An employee was deleted from the system.

EMP_EDIT -- Employee information was edited.

The report generated for each object instance will include the events that were defined by the System Administrator during the **Tracking Setup** function.


REPORT VIEW WINDOW

When a report has finished processing, it will appear in the **Report View** window. The user can scroll through the report, move from page to page, and zoom in on the report. In addition, the user can save the report to PDF, RTF, or text files, or print the report.

RAVEN-SB ASSESSMENT FORMS

The user can access the **Print** function for an individual assessment from the **RAVEN-SB Management** screen by clicking on the <PRINT> button located on the bottom of the screen. The Adobe Acrobat Reader will be launched with a PDF representation of the SB-MDS forms and a data overlay containing the assessment data. As a result, the *Print* dialog box will appear. Click on the <OK> button to print the complete assessment.

RAVEN-SB also has the capability to print multiple assessments contained in a standard export file. This capability can be accessed by selecting the Print Assessment(s) From File option from the OPTIONS menu item or the <PRINTER> icon located on the tool bar. The *RAVEN-SB Print Assessments* dialog box will be displayed.

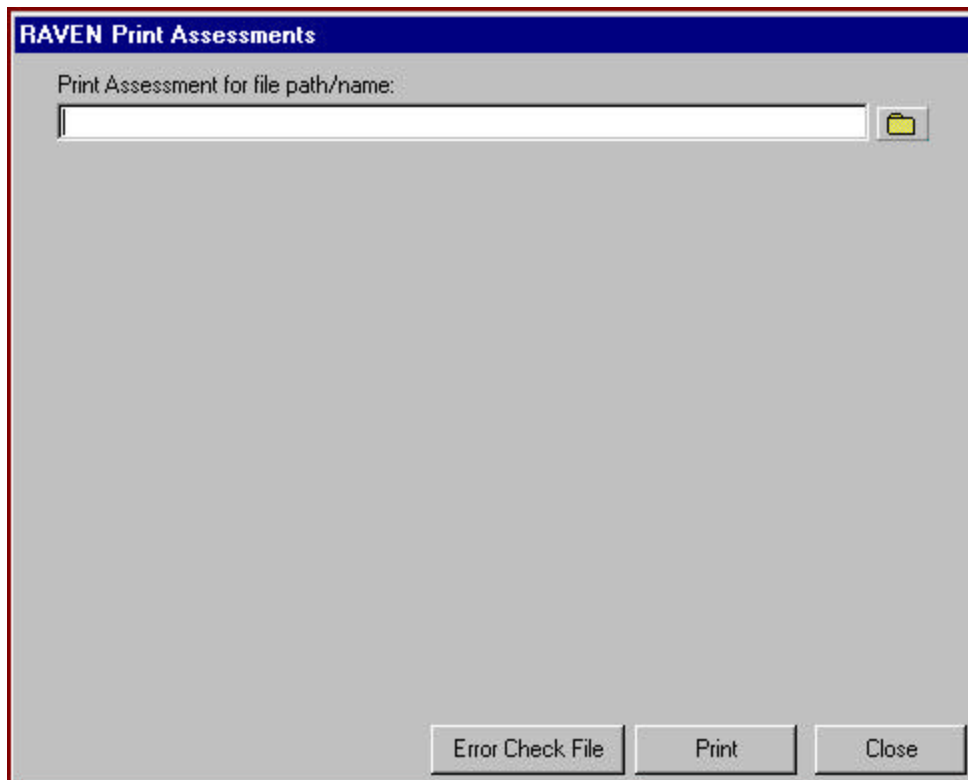

The user must provide the file path/name for the file containing the assessment(s). The user may then check for errors in the assessments by clicking on the <ERROR CHECK FILE> button. In order to proceed with printing, click on the <PRINT> button. A printer dialog box will appear for each assessment processed in the file. Click on the <OK> button to proceed. The user can stop printing assessments by pressing the <ESC> key. A message box will be displayed stating the number of assessments printed.

APPENDIX A

REPAIR/COMPACT UTILITY FOR MICROSOFT ACCESS DATABASE

REPAIR AND COMPACT


This utility module is named REPAIR.EXE. Use standard Windows procedures to execute this utility.



IMPORTANT: It is strongly recommended that the databases be backed up and archived before the Repair/Compact utility begins. Not all damaged Microsoft Access databases can be repaired and compacted by this or any repair/compact utility. If a database is badly corrupted, data truncation may occur as the result of a repair and compact and data loss may result.




IMPORTANT: This function only works on a valid Microsoft Access database with an ".MDB" or ".ARK" extension.

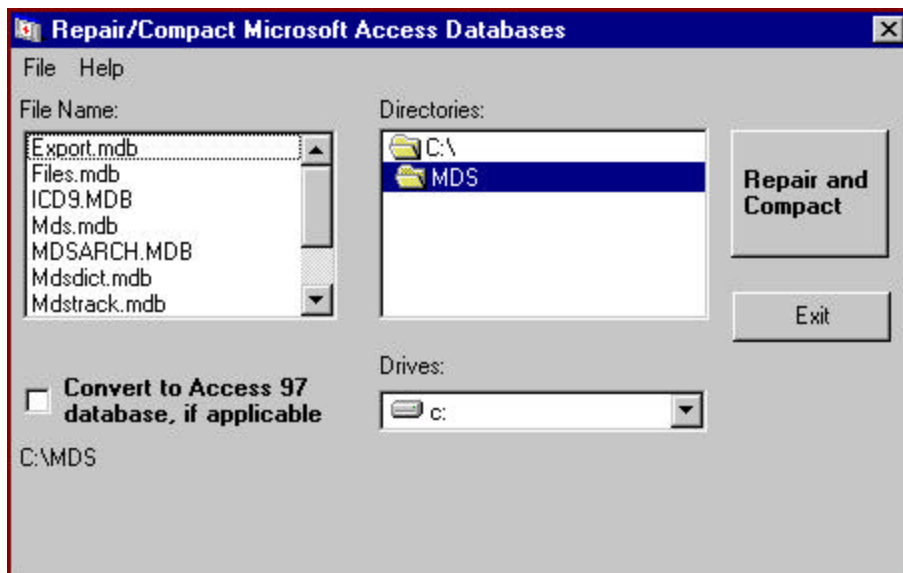

The Microsoft Access Databases Repair Utility is designed to fix a Microsoft Access database that becomes corrupted. To use the **Repair/Compact** utility, select the file using the drive, directory, and file list boxes on the single window that comprises this utility. The name of the file that has been selected appears in the lower left portion of the window.

The **Repair** function of the **Repair/Compact** utility tries to repair any index corruption and attempts to recover from data loss.

The **Compact** function of the **Repair/Compact** utility program compacts records on the selected database to recover space when records are deleted. Since the size of the database files does **not** decrease when records are deleted from the database, it is recommended that this program function be used when a large number of records have been deleted from the database.

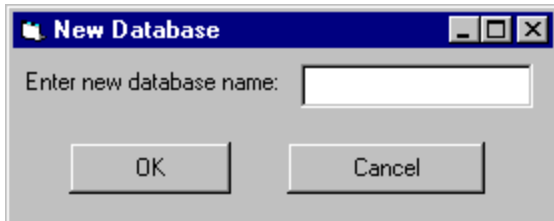

To begin the compacting process, click on the <REPAIR AND COMPACT> button. The user will be prompted for a new database name to which the data will be repaired and compacted. Enter the database name and click on the <OK> button to complete the **Repair/Compact** function. If successful, the message "Repair and compact completed!" is displayed otherwise an error message is displayed.

The amount of time it takes to repair and compact depends on the size of the database and the speed of the computer. To cancel the **Repair/Compact** function, click on the <CANCEL> button.



IMPORTANT: Do not interrupt the repairing and compacting process. An interruption may leave the database in an unstable state and could result in data loss.

After repairing, compacting, and verifying that the records have been repaired and compacted properly, it is important that the user replace the working database with the repaired and compacted one since the replacement is not performed automatically.

To exit the **Repair/Compact** utility, click on the <EXIT> button.

GLOSSARY OF TERMS

Data Entry

A MedQuest function that allows users to collect data for the data entry system designed by the **MedQuest Design** function by creating the data store file called **XXX.MDB** where "XXX" represents the module acronym (e.g., SWB.MDB is the data store file for the data entry module RAVEN-SB).

Data Entry Rule

A rule specified during the design process that is executed during data entry.

Data Entry Screen

A system designed by MedQuest to collect clinical data for a data analysis project. Each project is called a module and is represented by a three-character acronym (e.g., SWB is the Data Entry System for the Resident Assessment Validation and Entry for Swing Beds (RAVEN-SB) System project).

Module

See Data Entry System.

Screen

A screen is an area beneath the tab where variables are presented. A tab represents a screen. The variables on the screen are displayed by selecting a tab. A tab may contain one or more subtabs.

Subscreen

A subscreen is an area beneath the subtab where variables are presented. A subtab represents a subscreen. The variables on the subscreen are displayed by selecting a subtab.

Tab/Subtab

A tab/subtab indicates the label for the screen that can be selected (e.g., Tab **Section A** is the tab used for retrieving the **Section A** screen).

Variable

A variable is a data entry field that accepts data according to the specifications and rules indicated by the designer during the design.

Variable Type

Variable types and attributes determine the type of data the system will accept (e.g., a date variable will accept only data in date format, i.e., MM/DD/YYYY) during data entry. The designer defines a variable type for each variable. For example, if variable **Admission Date** accepts only data that are in date form, the variable type must be defined as variable type.